


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

4-Benzoyloxy-2-azetidinone

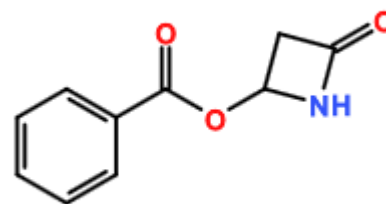

## CHROMATOGRAM

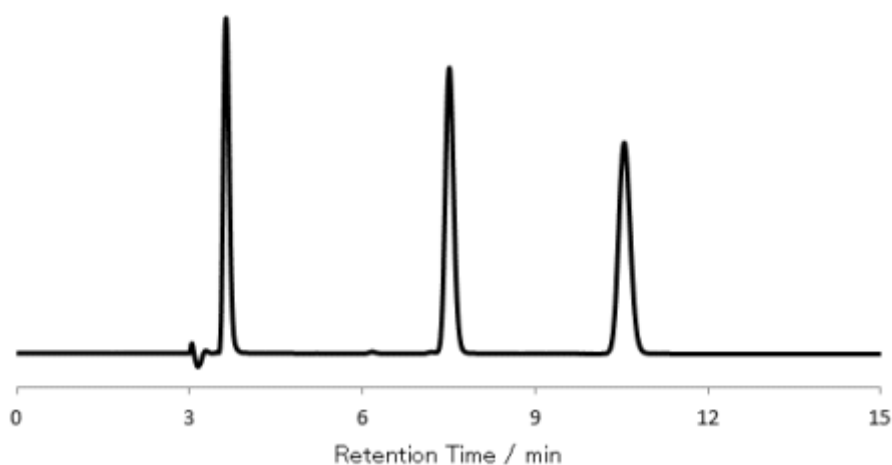

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IK		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / 2-propanol = 60 / 40		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	5.0	Sample Concentration	1.0 mg/mL in EtOH
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.5	1.50	1 - 2	9.23
2	10.5	2.51		