


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Ethiazide

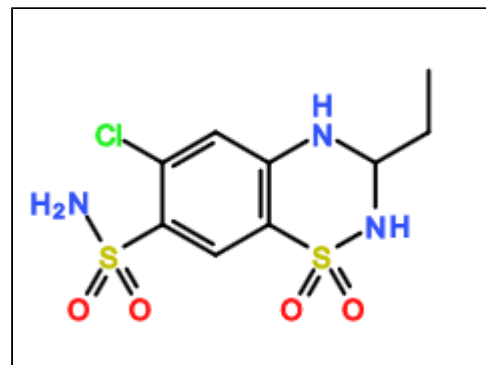

## CHROMATOGRAM

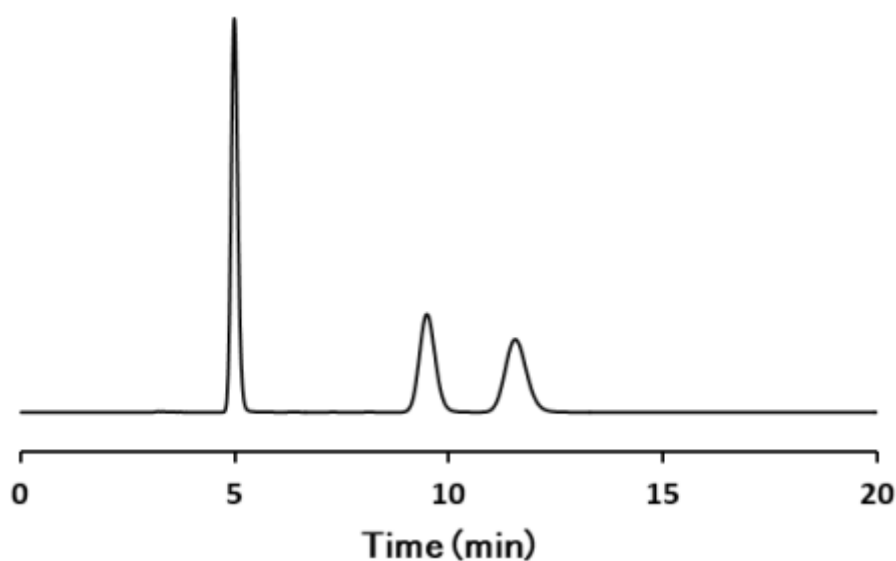

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IK		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	100mM KPF6 aq. / acetonitrile = 80 / 20		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	5	Sample Concentration	1.0 mg/mL
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.5	2.16	1 - 2	2.53
2	11.6	2.85		