


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

N-CBZ-DL-Valine

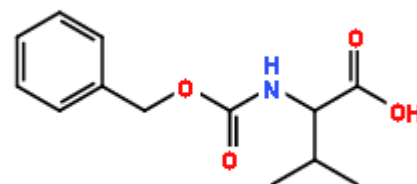

CHROMATOGRAM

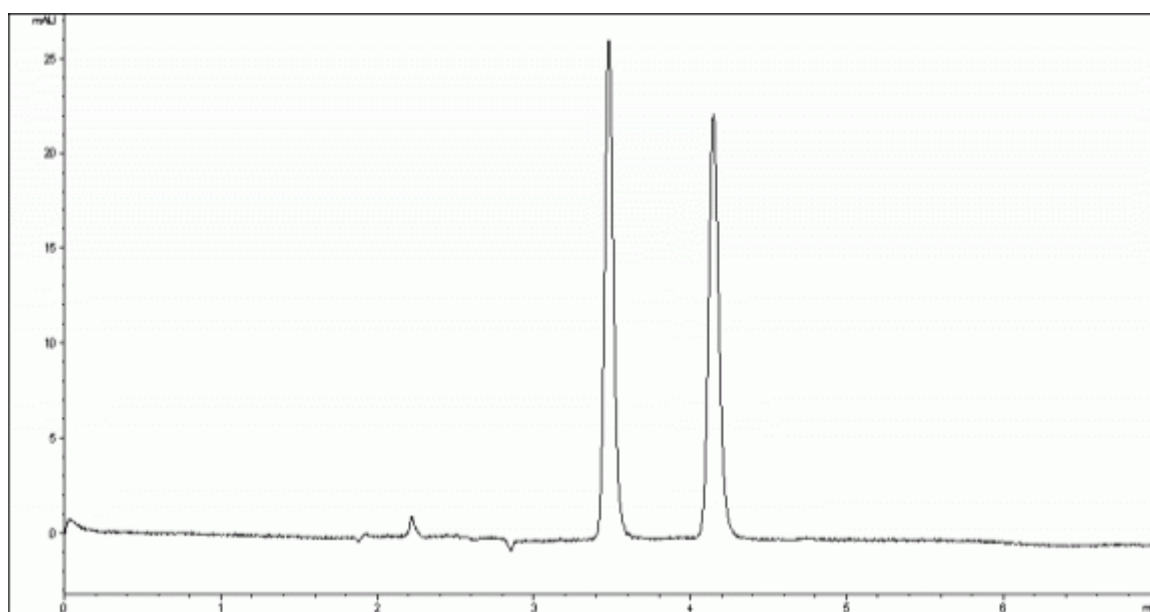

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	n-hexane / chloroform / ethanol / trifluoroacetic acid = 68 / 30 / 2 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.5	0.92	1 - 2	6.10
2	4.2	1.28		