


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Pranoprofen

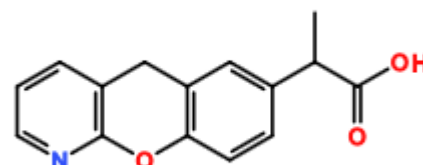

## CHROMATOGRAM

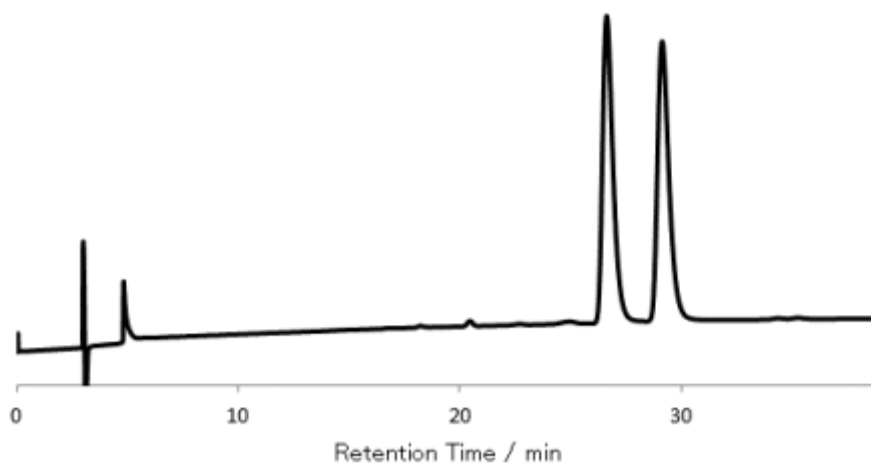

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IK-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / 2-propanol / trifluoroacetic acid = 98 / 2 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	5	Sample Concentration	1 mg/mL in Mobile phase
Detection	UV-VIS 248 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	26.6	7.88	1 - 2	1.11
2	29.1	8.71		2.83
3				