


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Tulobuterol

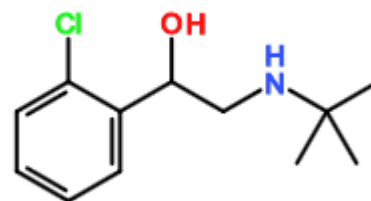

CHROMATOGRAM

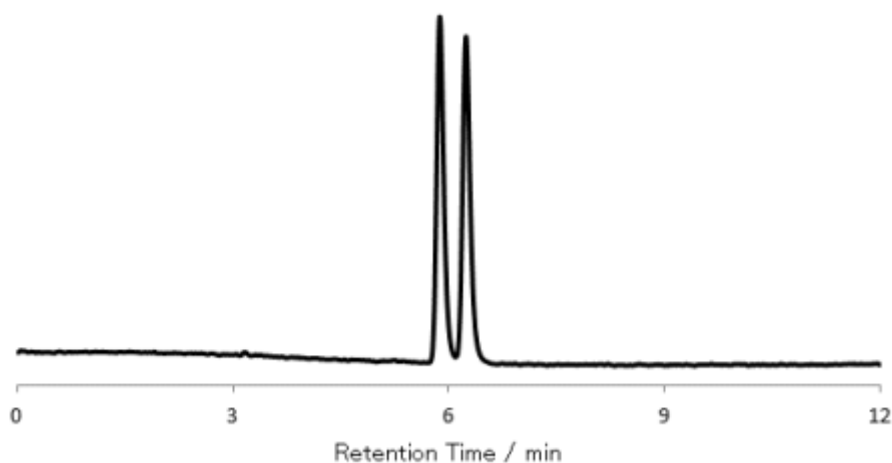

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IM		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / 2-propanol / diethylamine = 98 / 2 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	10	Sample Concentration	0.5 mg/mL in Mobile phase
Detection	UV-VIS 263 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.9	0.96	1 - 2	1.76
2	6.2	1.08		