


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Ornidazole

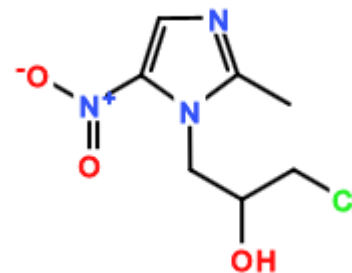

## CHROMATOGRAM

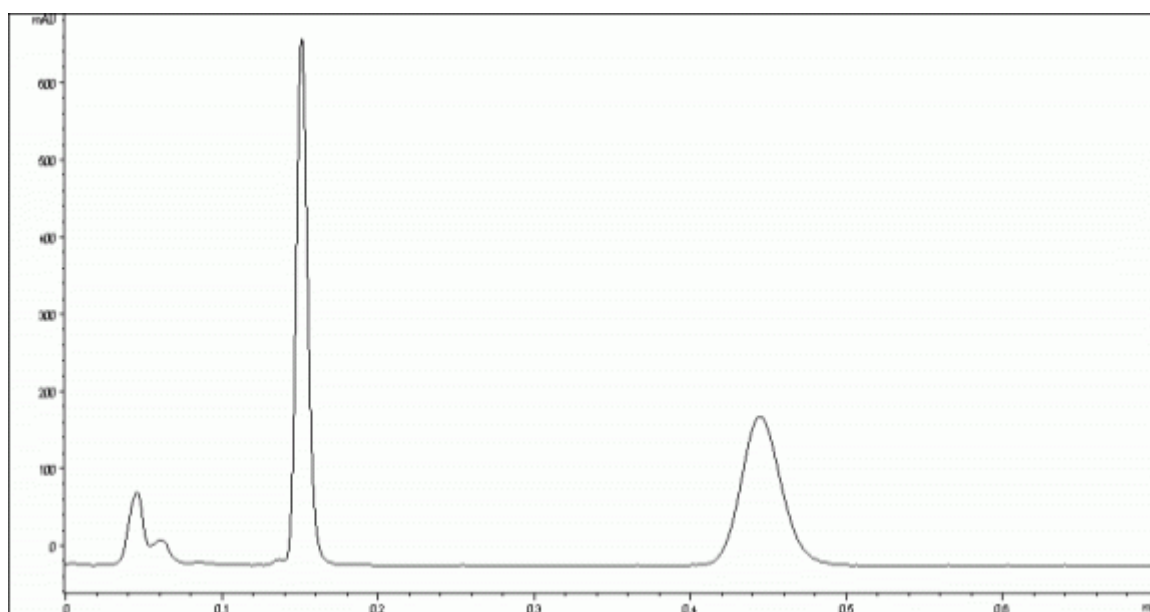

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 50 mm (3 µm)		
Mobile Phase Composition	dichloromethane / methanol / ethanolamine = 88 / 12 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	5	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	0.2	0.31	1 - 2	9.64
2	0.5	2.87		