


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Oxprenolol

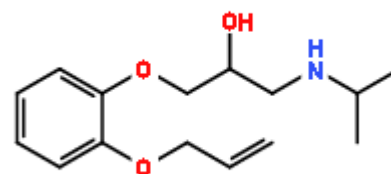

CHROMATOGRAM

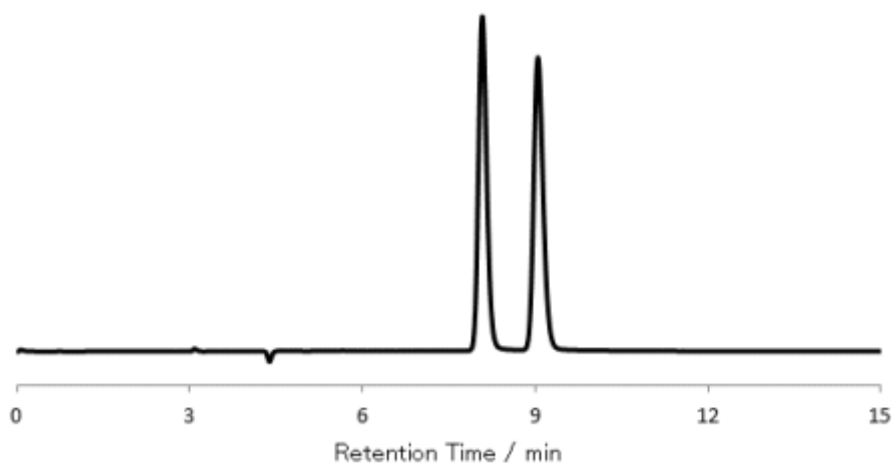

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IN		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 93 / 7 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	5	Sample Concentration	1.0 mg/mL in Mobile phase
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	8.1	1.69	1 - 2 1.19	3.27
2	9.0	2.01		