


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

p-Anisoic

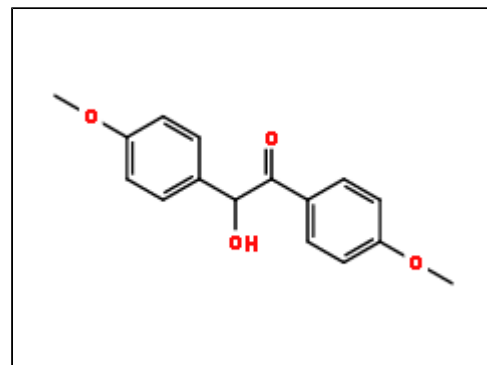

CHROMATOGRAM

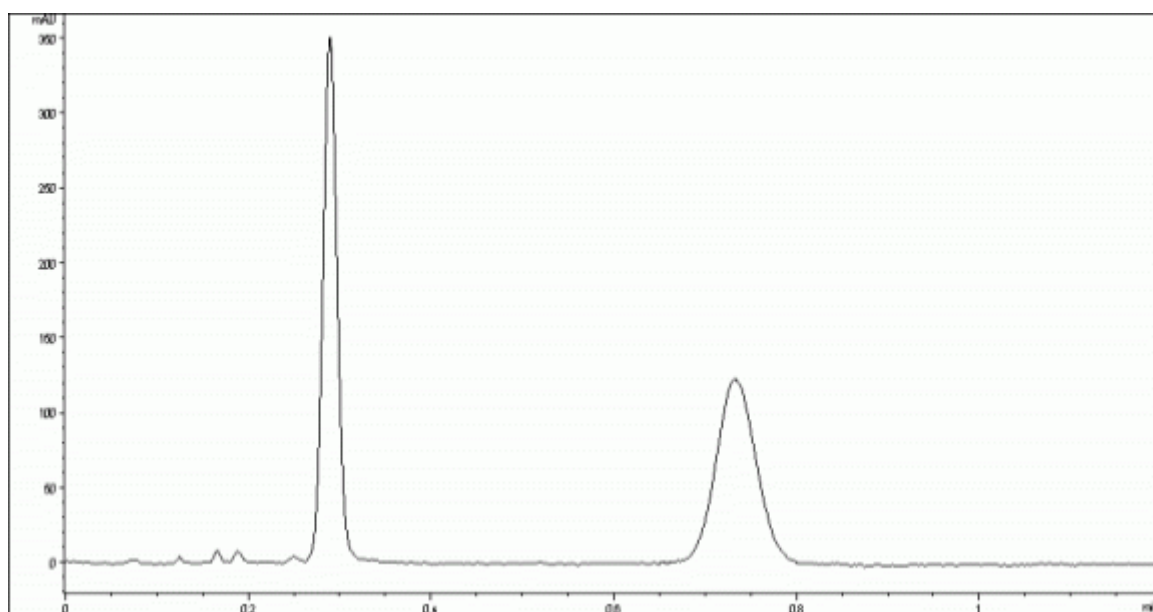

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 50 mm (3 µm)		
Mobile Phase Composition	methyl tert-butyl ether / ethanol = 90 / 10		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	5	Temperature (°C)	25
Injection Volume (µL)	0.002 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	0.3	1.52	1 - 2	8.07
2	0.7	5.37		