


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Indapamide

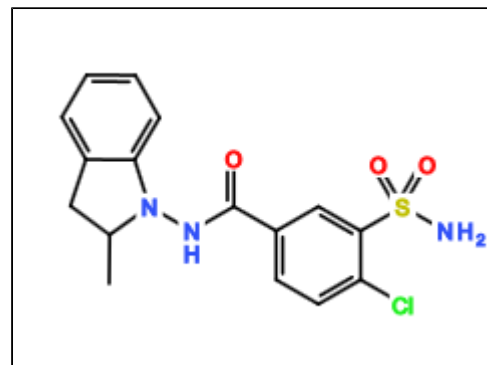

CHROMATOGRAM

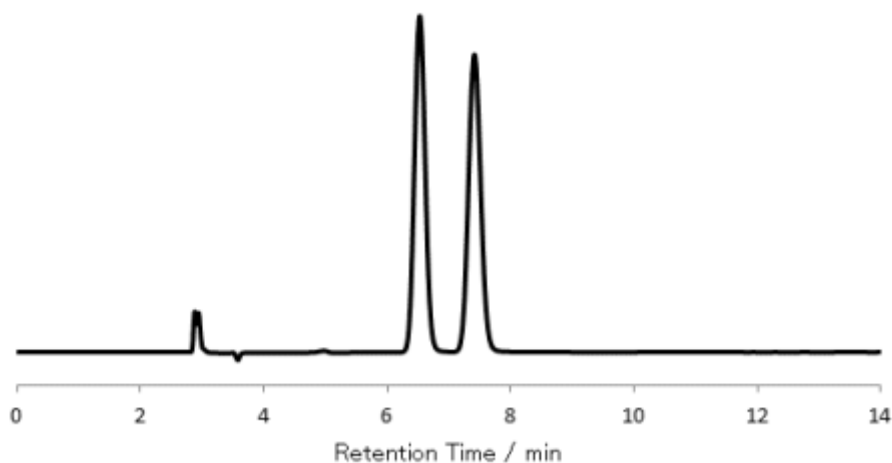

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IN		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / tetrahydrofuran / diethylamine = 50 / 50 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	10	Sample Concentration	0.5 mg/mL
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.5	1.17	1 - 2	2.49
2	7.4	1.47		