


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Indapamide

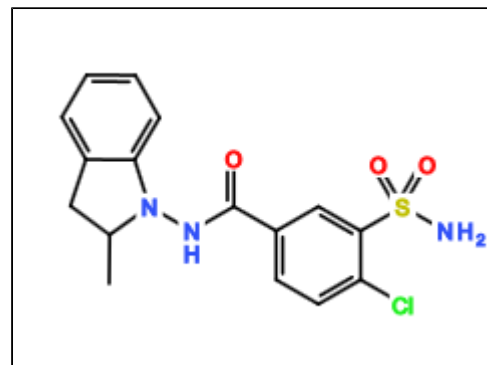

## CHROMATOGRAM

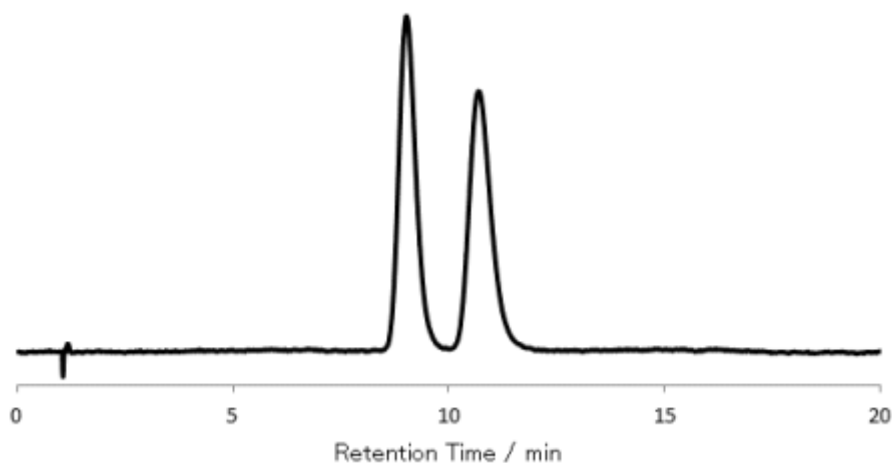

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IN		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	CO <sub>2</sub> / methanol / isopropylamine = 80 / 20 / 0.02		
Chromatographic Mode	SFC		
Flow rate (mL/min)	3.0	Temperature (°C)	40
Injection Volume (µL)	5	Sample Concentration	1.0 mg/mL in Mobile phase
Detection	UV-VIS 240 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.0	7.20	1 - 2	2.03
2	10.7	8.72		