


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Pemetrexed disodium

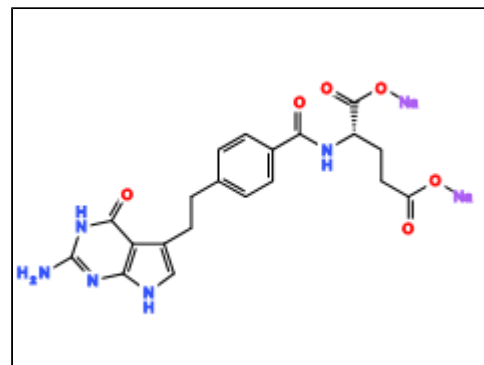

## CHROMATOGRAM

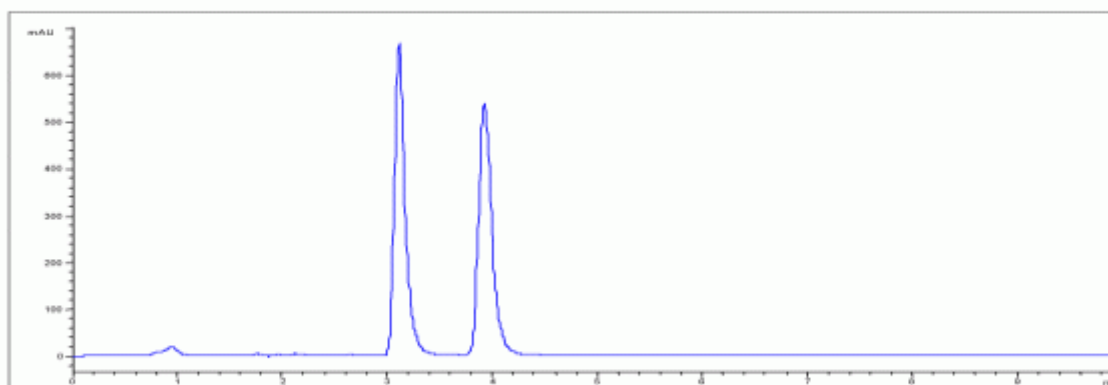

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	dichloromethane / methanol / trifluoroacetic acid = 90 / 10 / 0.5		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 240 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.2	0.78	1 - 2	4.03
2	4.1	1.29		