


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

(±)-N,N'-Bis(alpha-methylbenzyl) sulfamide

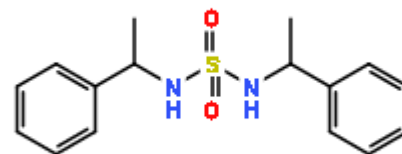

## CHROMATOGRAM

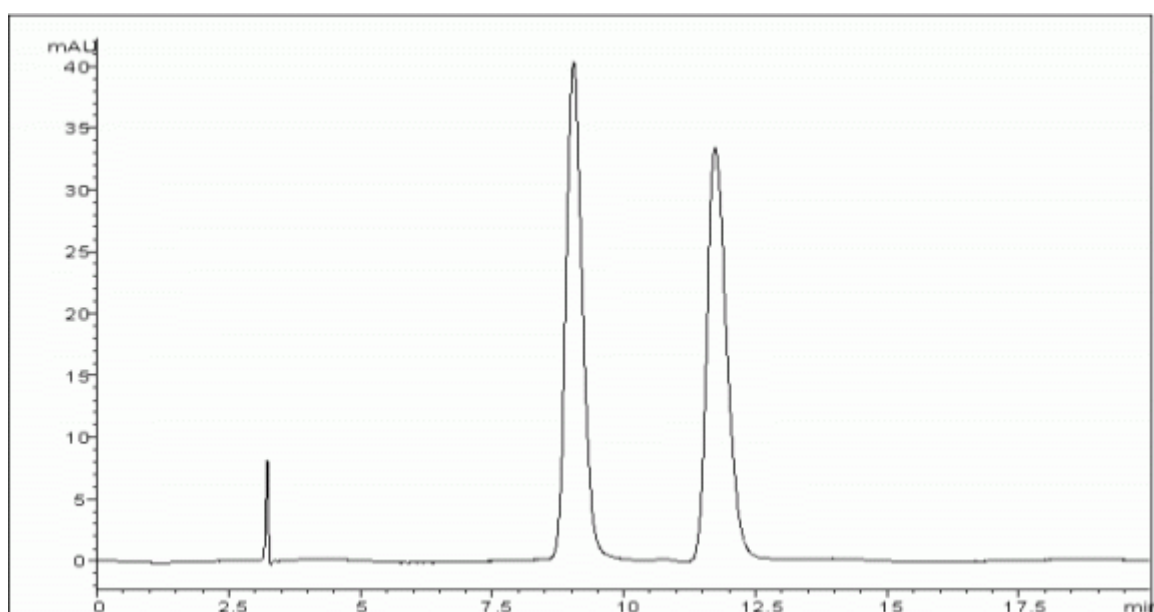

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	methyl tert-butyl ether = 100		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.0	1.95	1 - 2	4.36
2	11.7	2.83		