


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

2-Ethyl-1-Indanone

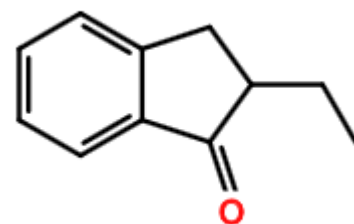

## CHROMATOGRAM

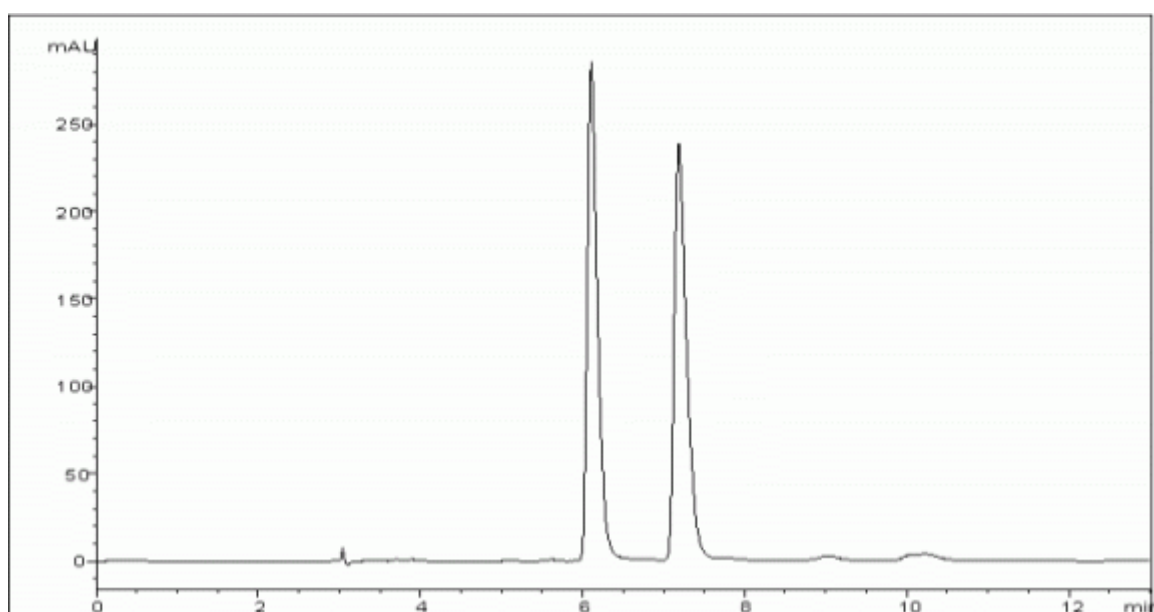

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / methyl tert-butyl ether = 50 / 50		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.1	1.09	1 - 2	4.38
2	7.2	1.46		