


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

2-Phenoxypropionic acid

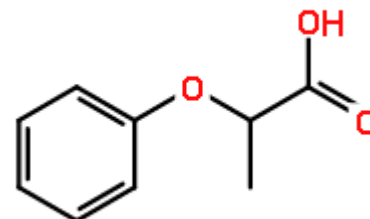

CHROMATOGRAM

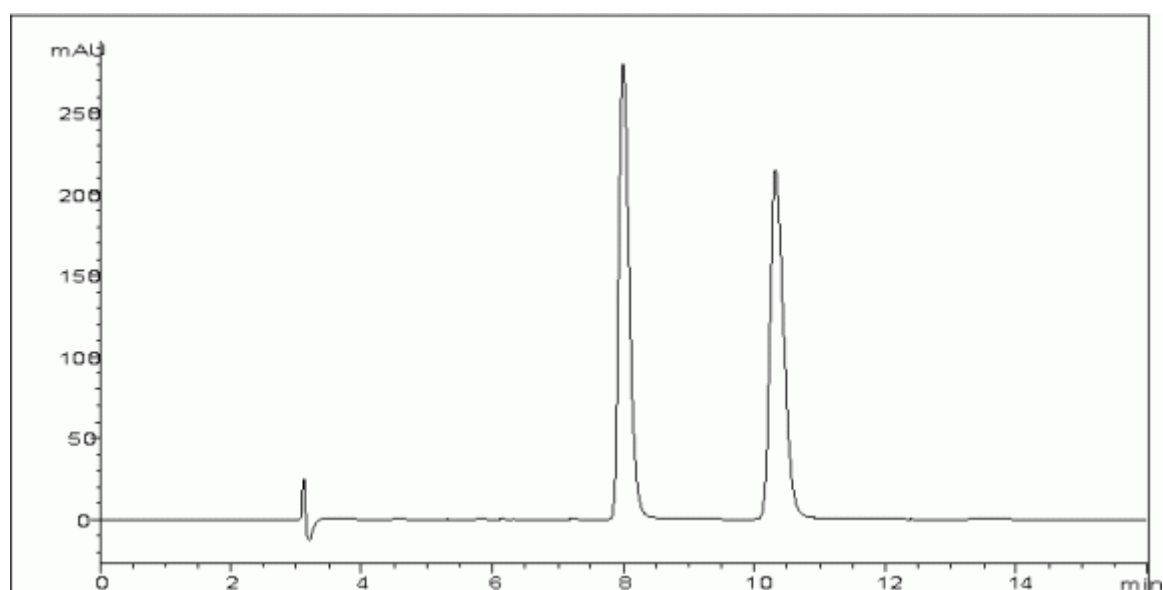

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / dimethoxymethane / trifluoroacetic acid = 75 / 25 / 0.1		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	8.0	1.68	1 - 2	7.03
2	10.3	2.46		