


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

alpha-(2,4-Dichlorophenyl)-1H-imidazole-1-ethanol

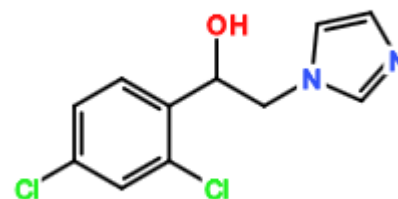

CHROMATOGRAM

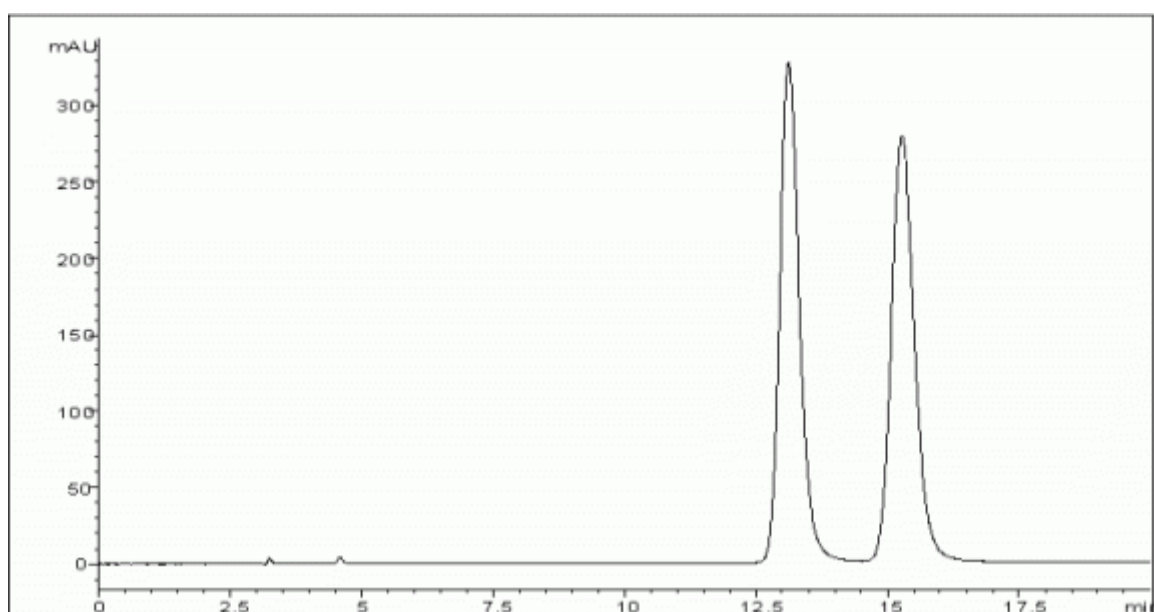

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-heptane / 2-propanol / diethylamine = 85 / 15 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	13.1	3.37	1 - 2	2.97
2	15.3	4.09		