


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Amlodipine

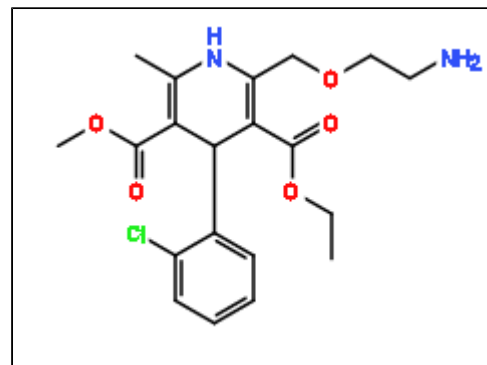

CHROMATOGRAM

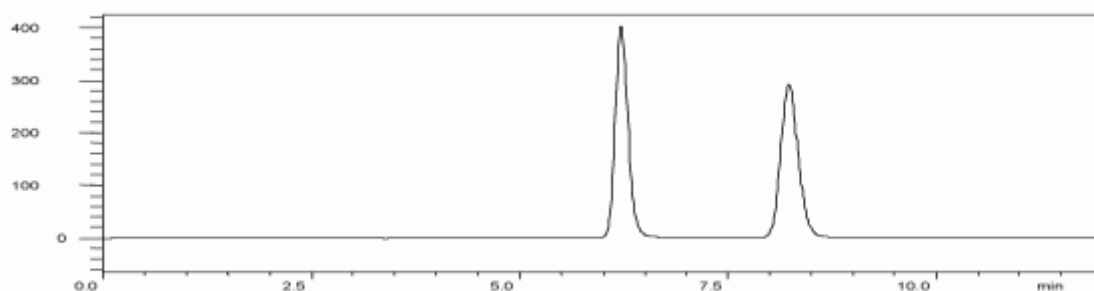

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	dichloromethane / acetonitrile / isopropylamine = 85 / 15 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	15
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 360 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.2	1.07	1 - 2	5.52
2	8.2	1.74		