


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Atomoxetine (racemic)

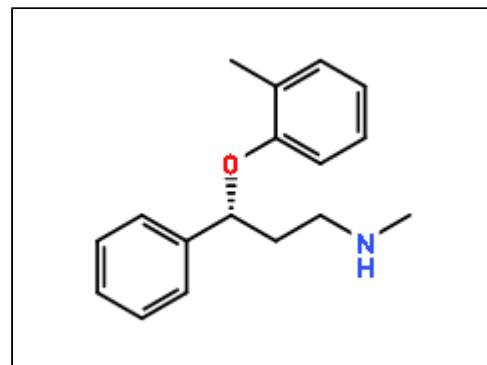

CHROMATOGRAM

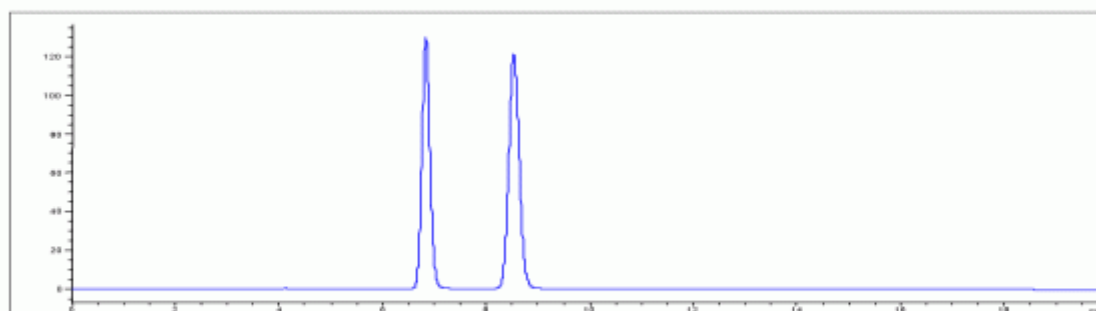

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / 2-propanol / trifluoroacetic acid / diethylamine = 80 / 20 / 0.1 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.5	1.17	1 - 2	5.02
2	8.5	1.84		