


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

1-Methyl-DL-tryptophan

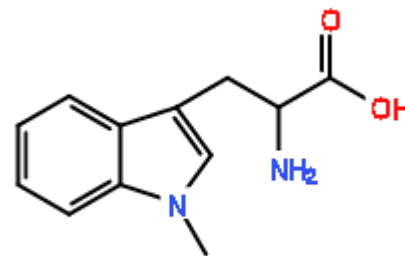

CHROMATOGRAM

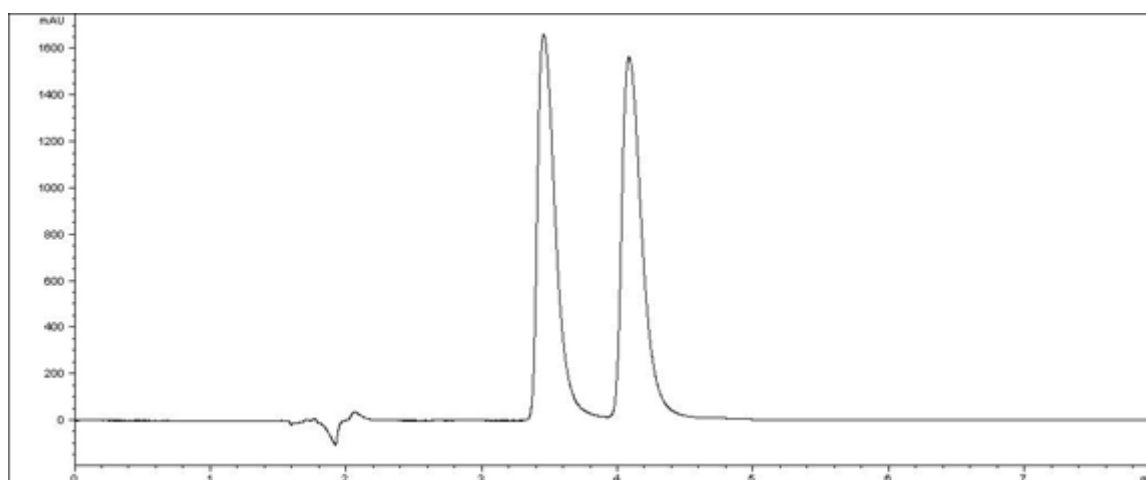

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® ZWIX(+) 3 x 150 mm (3 µm)		
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / water = 90 / 10		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order	D/L		

No.	Rt	K'	Alpha Value	Resolution
1	3.5	1.08	1 - 2	2.53
2	4.1	1.46		