


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Cetirizine

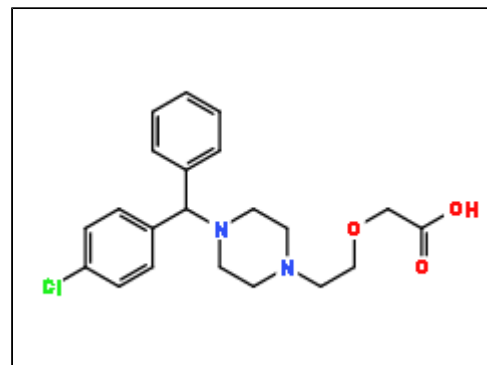

## CHROMATOGRAM

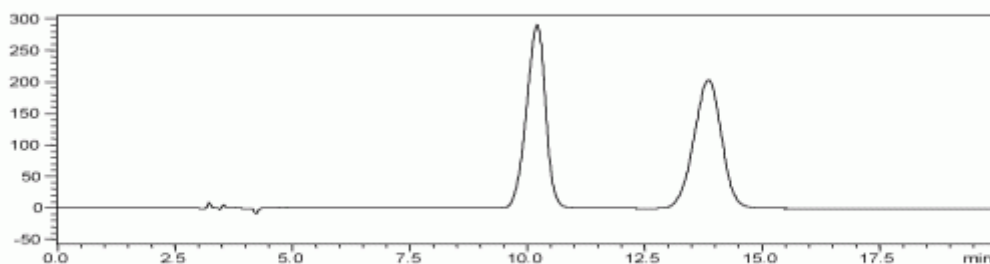

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / 1,4-dioxane / ethylenediamine = 40 / 50 / 10 / 0.3		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	30
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	10.2	2.40	1 - 2	3.90
2	13.8	3.60		