


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Diperodon

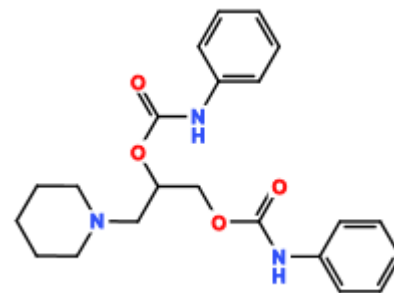

CHROMATOGRAM

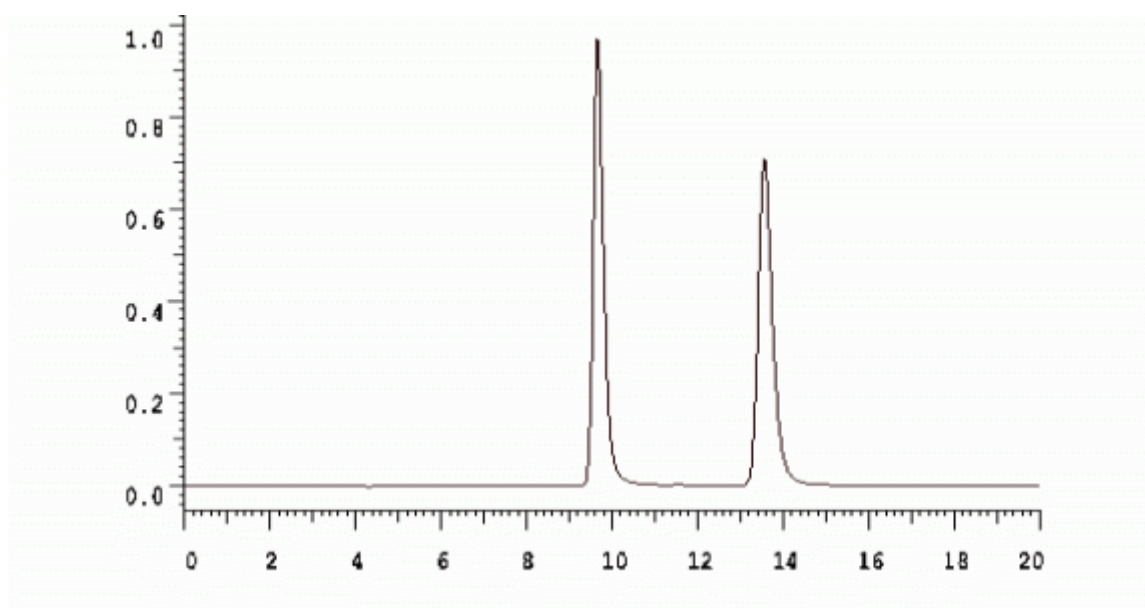

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / dichloromethane / diethylamine = 50 / 50 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 240 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.7	2.28	1 - 2	7.68
2	13.6	3.61		