


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Fenoxaprop-P-ethyl

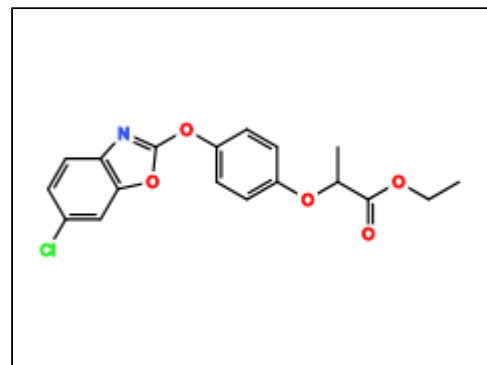

CHROMATOGRAM

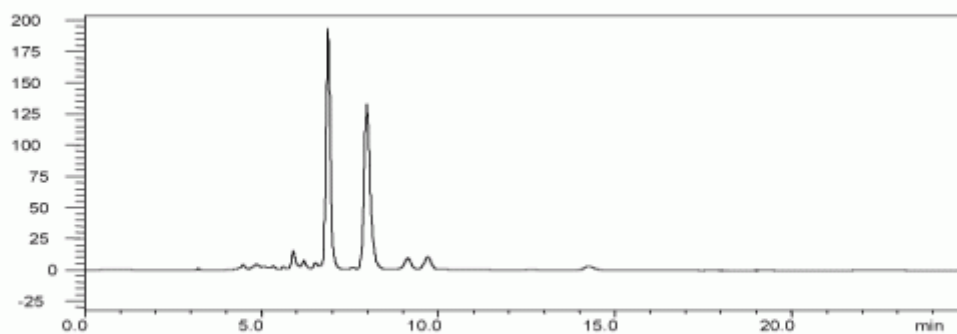

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 80 / 20 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 235 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.8	1.27	1 - 2	3.82
2	7.9	1.64		