


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Laudanosine

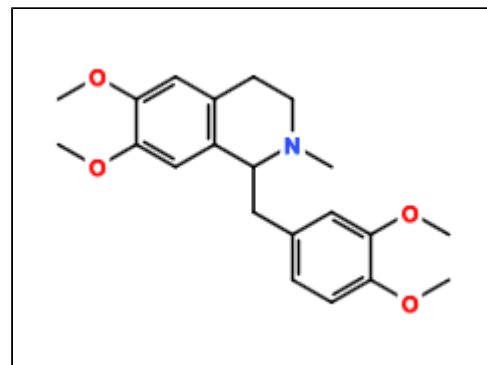

## CHROMATOGRAM

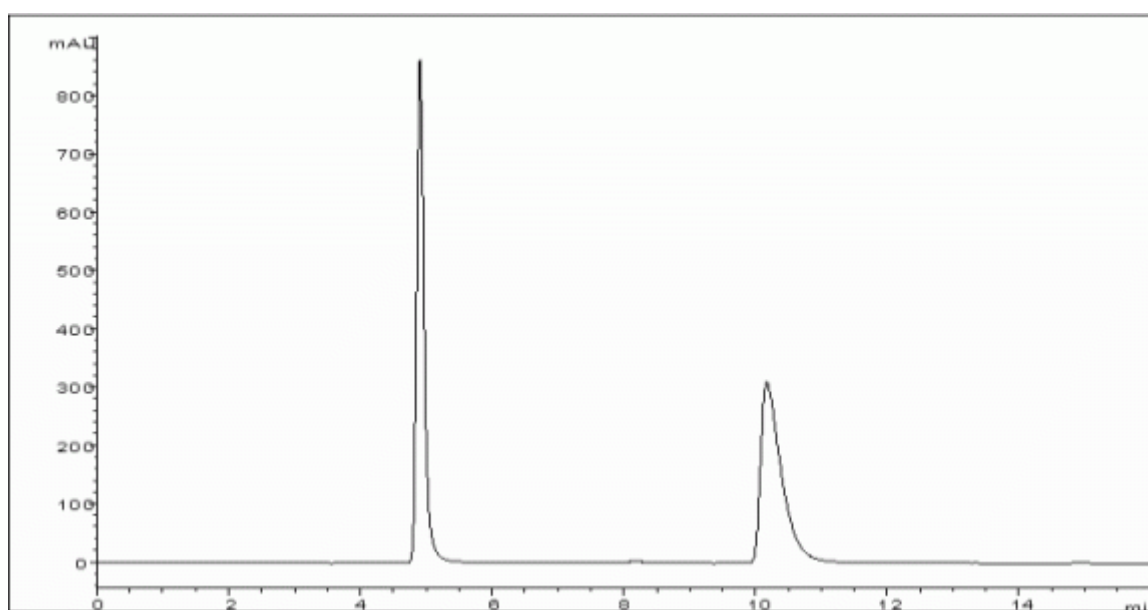

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	methyl tert-butyl ether / ethanol / ethylenediamine = 95 / 5 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 290 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.9	0.64	1 - 2	14.24
2	10.2	2.40		