


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Metanephrine

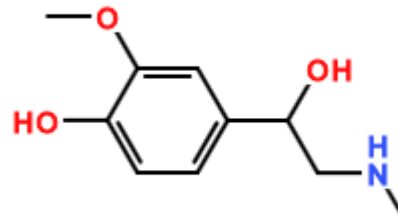

## CHROMATOGRAM

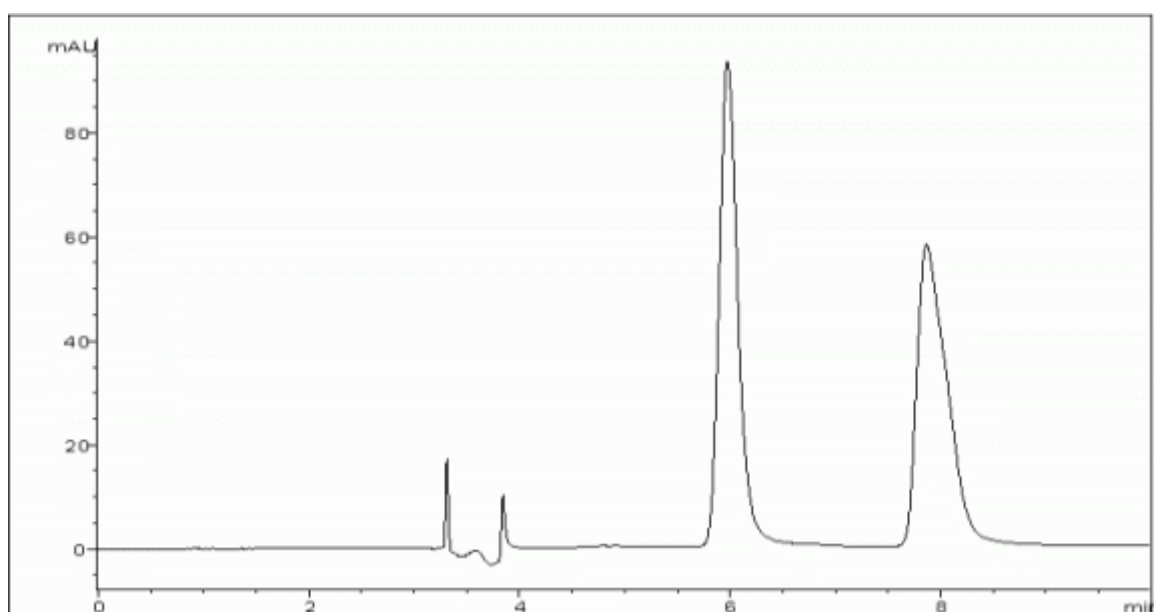

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	methyl tert-butyl ether / methanol / ethanolamine = 94 / 6 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.0	0.95	1 - 2	4.48
2	7.9	1.57		