


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Metanephrine

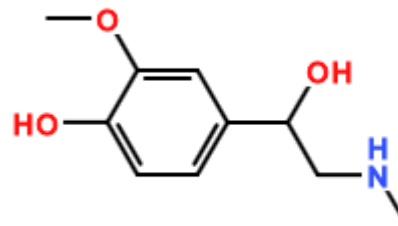

## CHROMATOGRAM

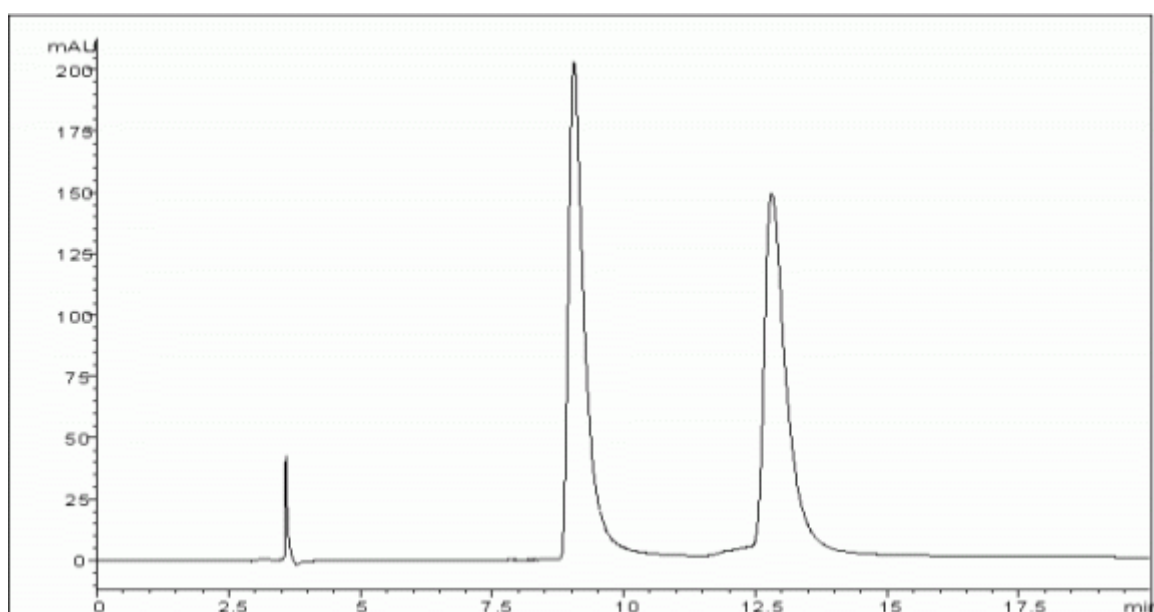

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / 2-propanol / diethylamine = 70 / 30 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.1	1.98	1 - 2	6.17
2	12.8	3.23		