


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

N-acetyl-DL-Phenylalanine

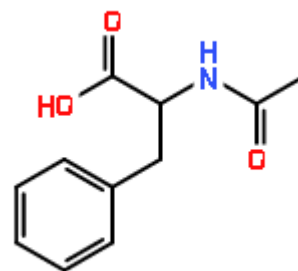

## CHROMATOGRAM

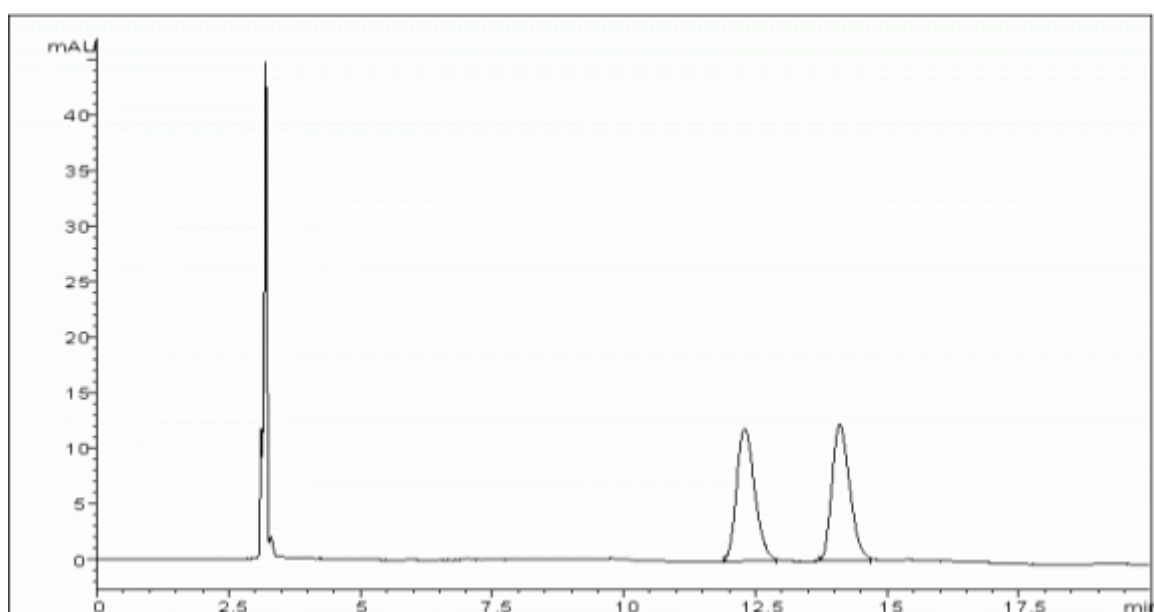

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / tetrahydrofuran / trifluoroacetic acid = 75 / 25 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	12.3	3.02	1 - 2	2.92
2	14.1	3.62		