


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

N-CBZ-DL-Alanine

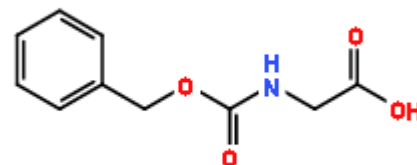

## CHROMATOGRAM

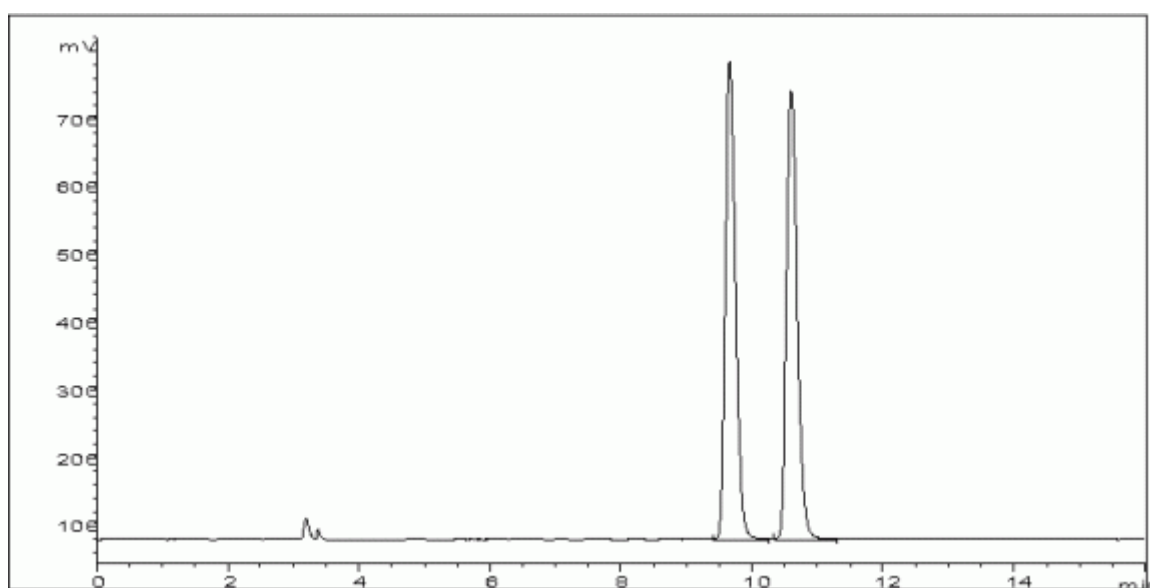

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethyl acetate / trifluoroacetic acid = 70 / 30 / 0.1		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.008 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.7	2.16	1 - 2	3.21
2	10.6	2.46		