


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

DL-Ala-DL-Leu-Gly

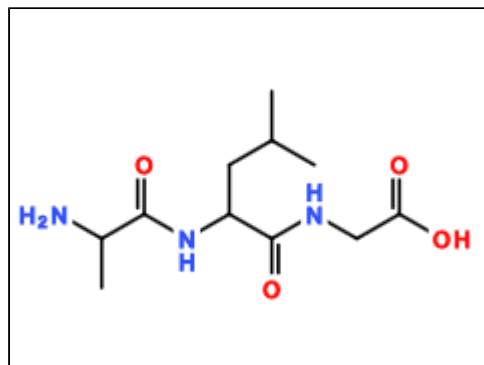

## CHROMATOGRAM

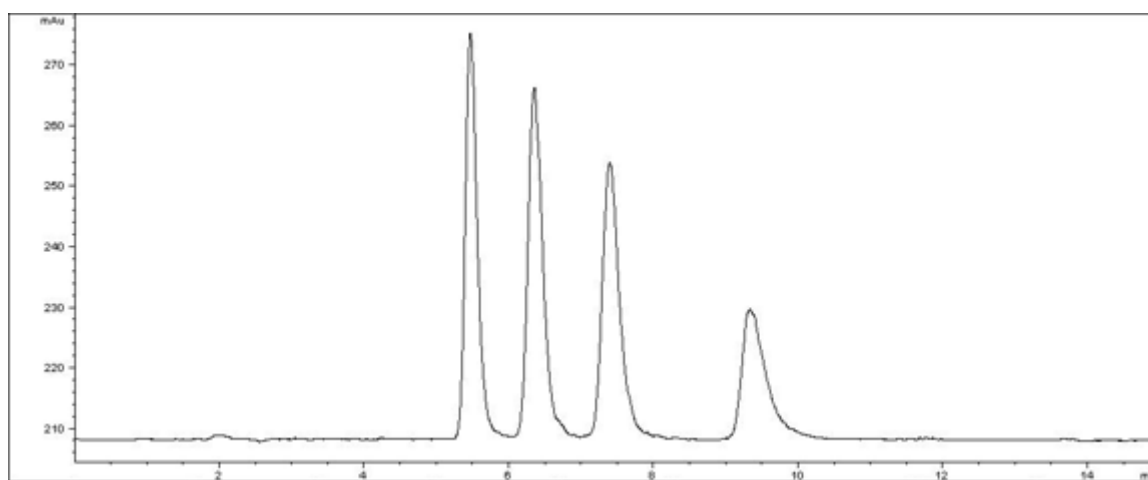

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(-)		
Column Size (Particle Size)	3 x 150 mm (3 µm)		
Mobile Phase Composition	50mM formic acid + 25mM diethylamine in methanol / tetrahydrofuran / water = 49 / 49 / 2		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.0015 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.5	2.08	1 - 2	2.71
2	6.4	2.58	2 - 3	2.66
3	7.4	3.17	3 - 4	3.92
4	9.3	4.26		