


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Paliperidone

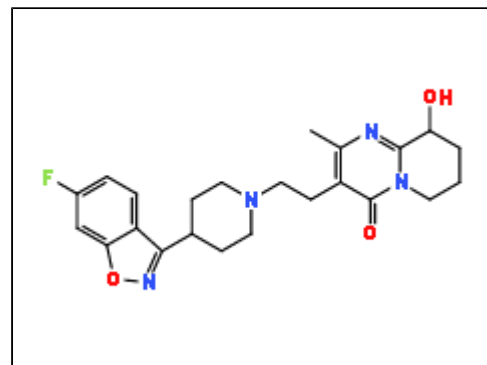

## CHROMATOGRAM

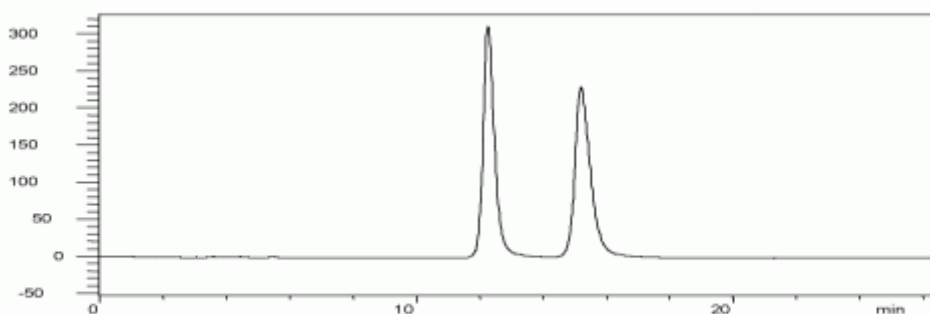

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 50 / 50 / 0.3		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 280 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	12.2	3.07	1 - 2	1.32
2	15.1	4.04		