


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Praziquantel

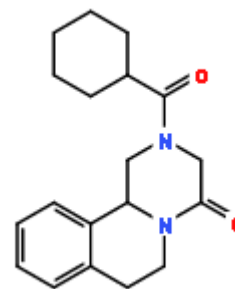

## CHROMATOGRAM

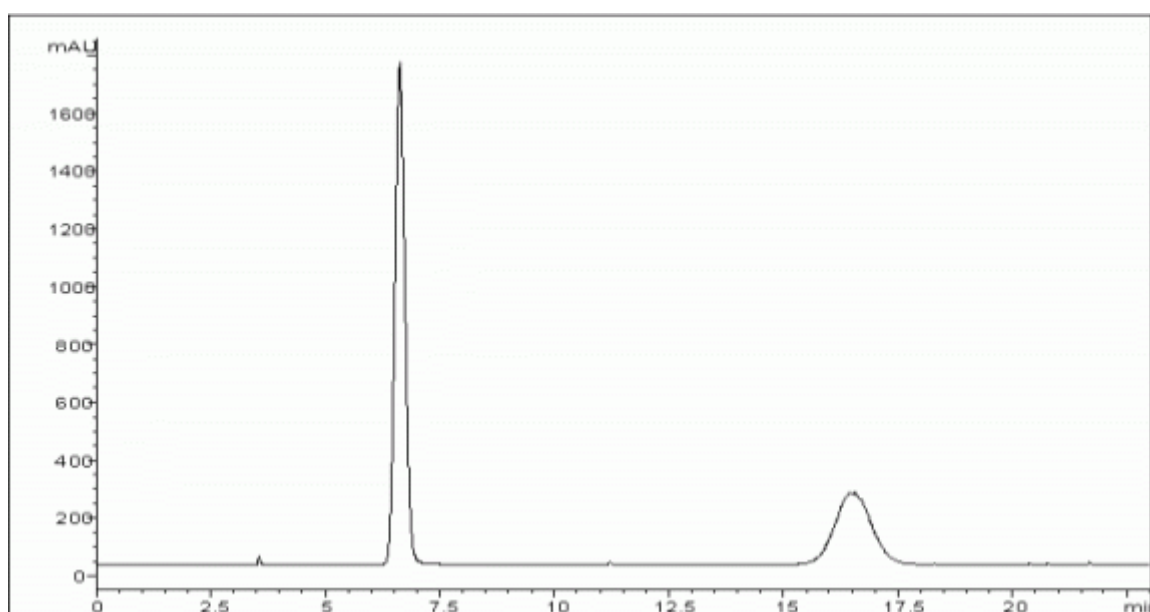

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	dichloromethane / methanol = 98 / 2		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	ELSD		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.6	1.16	1 - 2	3.76
2	16.5	4.36		10.49