


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Rabeprazole

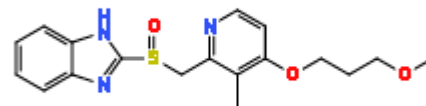

CHROMATOGRAM

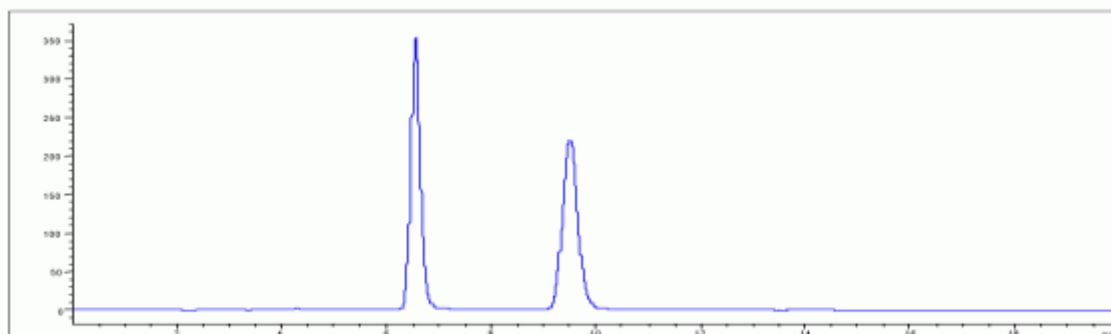

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 30 / 70 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	35
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 280 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.6	1.21	1 - 2	6.04
2	9.5	2.17		