


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Rosuvastatin (racemic)

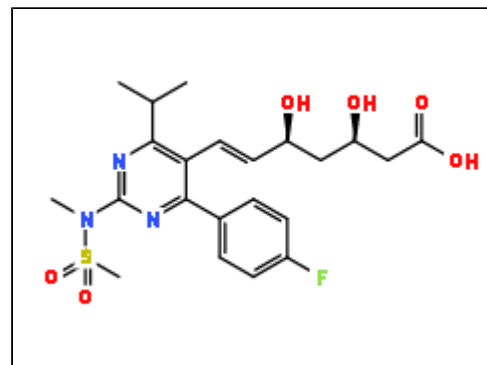

CHROMATOGRAM

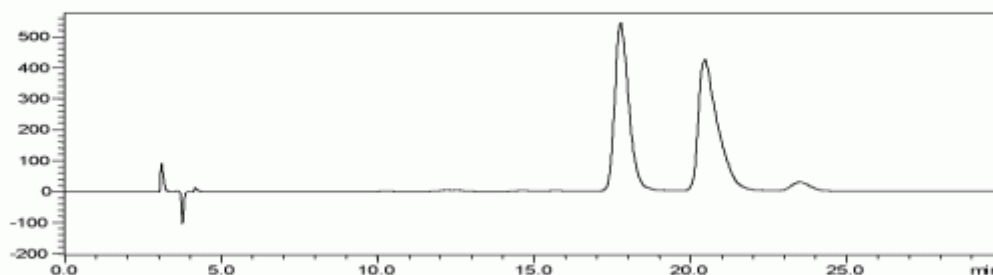

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / 2-propanol / trifluoroacetic acid = 90 / 5 / 5 / 0.3		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 250 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	17.7	4.91	1 - 2	2.30
2	20.4	5.81		