


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Solifenacin (containing two pairs of enantiomers)

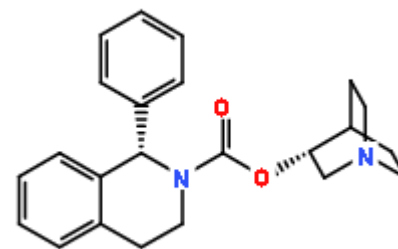

CHROMATOGRAM

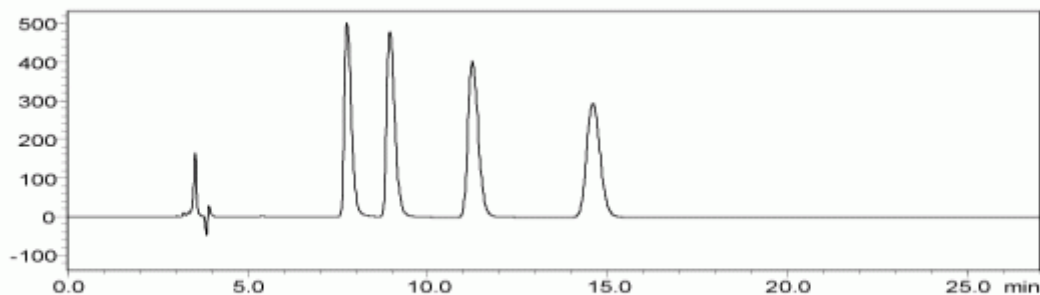

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / 2-propanol / ethanol / diethylamine = 60 / 25 / 15 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	30
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.8	1.60	1 - 2	1.27
2	9.1	2.03	2 - 3	1.40
3	11.5	2.84	3 - 4	1.39
4	14.8	3.94		