


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Tramadol

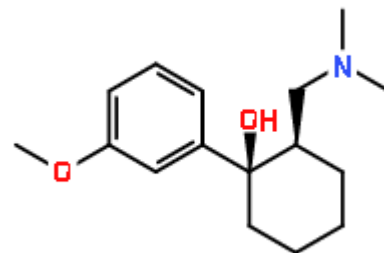

CHROMATOGRAM

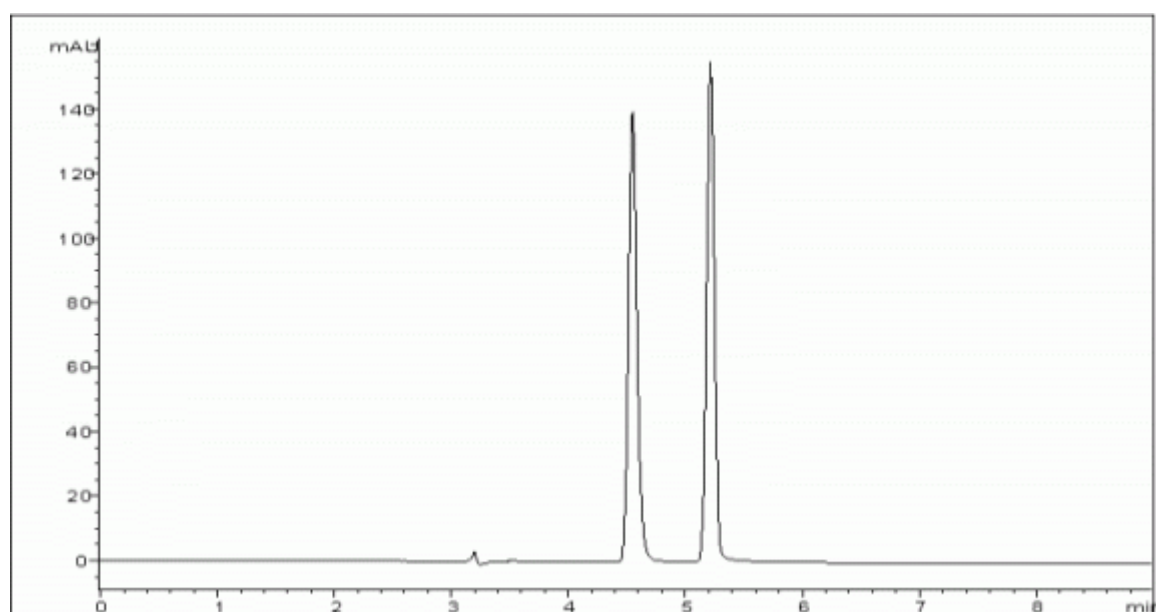

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / methyl tert-butyl ether / ethanol / ethanolamine = 69.5 / 29.5 / 1 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.003 (mg)	Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.6	0.48	1 - 2	5.23
2	5.2	0.70		