


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

trans-Stilbene Oxide

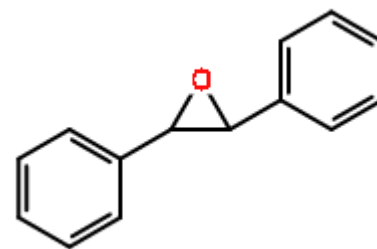

CHROMATOGRAM

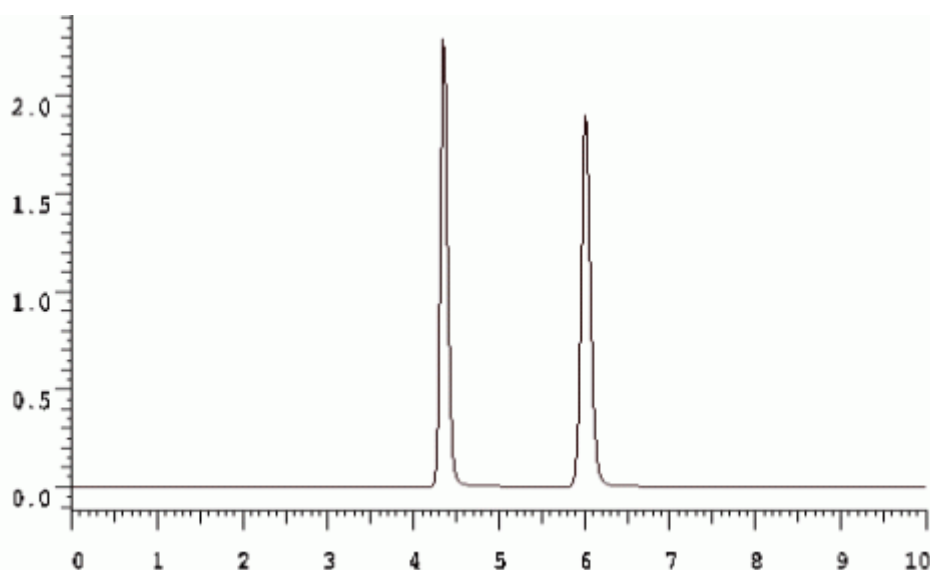

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / dichloromethane = 75 / 25		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 240 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.4	0.48	1 - 2	9.92
2	6.0	1.04		