


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Cyclandelate (mixture of 2 pairs of enantiomers)

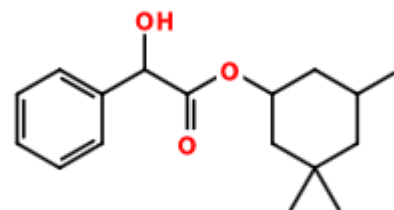

## CHROMATOGRAM

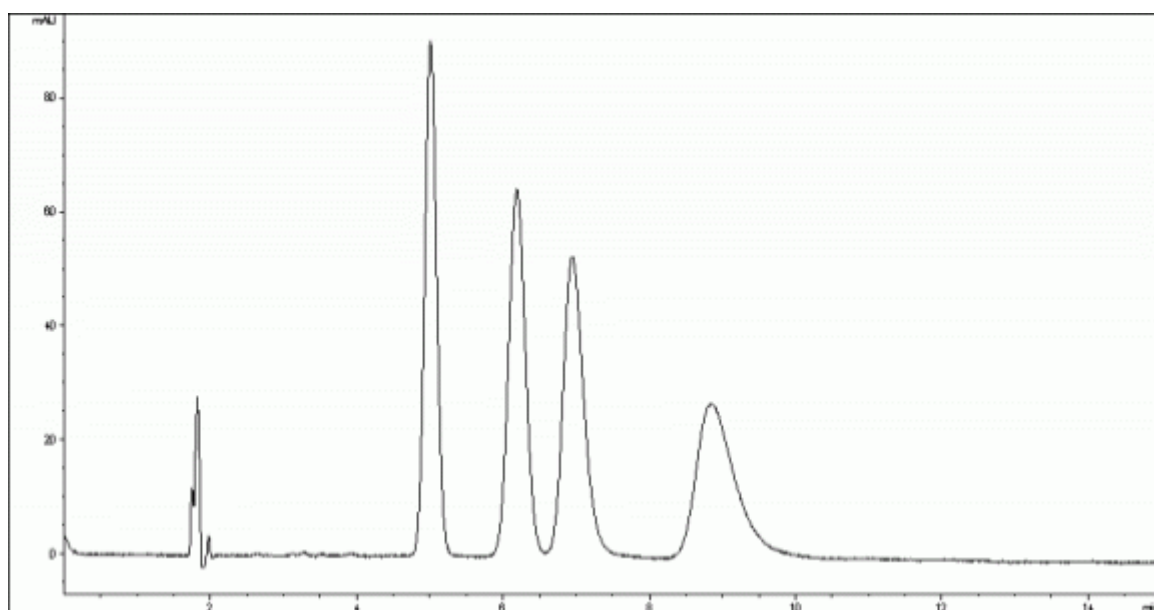

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AY-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol = 90 / 10		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.0	1.82	1 - 2	1.36
2	6.2	2.49	2 - 3	1.17
3	6.9	2.92	3 - 4	1.37
4	8.9	4.00		