


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Ethyl alpha-chlorophenylacetate

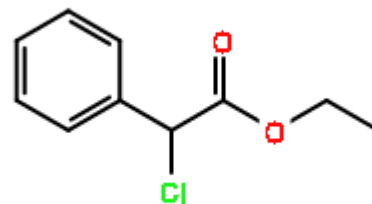

CHROMATOGRAM

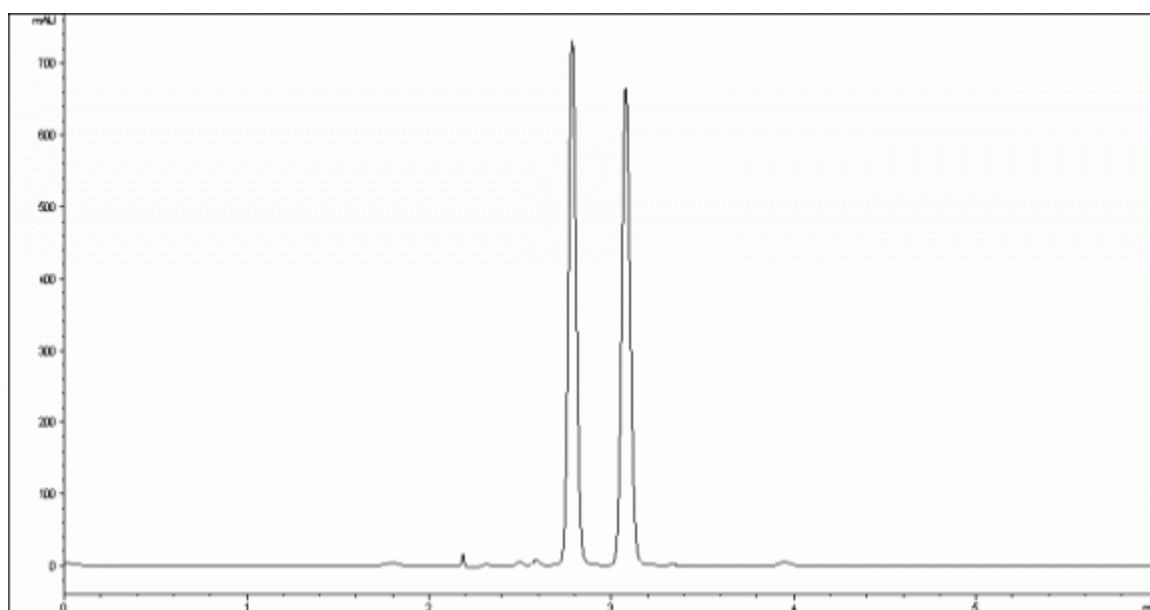

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AY-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	n-hexane / 2-propanol = 80 / 20		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.8	0.57	1 - 2	3.65
2	3.1	0.74		