


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Nomifensine

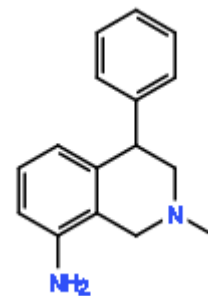

## CHROMATOGRAM

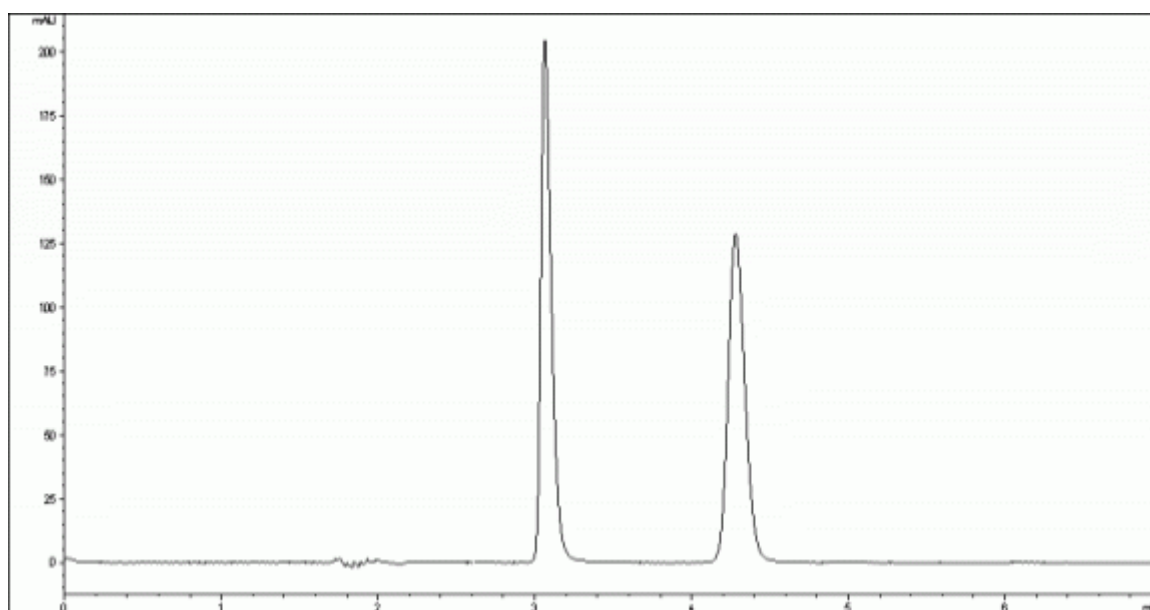

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK <sup>®</sup> AY-3		
Column Size (Particle Size)	4.6 x 150 mm (3 $\mu$ m)		
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 80 / 20 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature ( $^{\circ}$ C)	25
Injection Volume ( $\mu$ L)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.1	0.73	1 - 2	7.84
2	4.3	1.42		