


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Pindolol

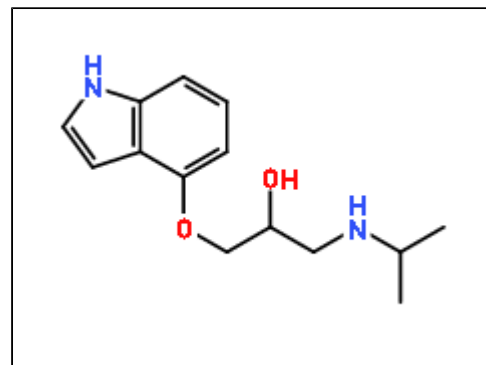

CHROMATOGRAM

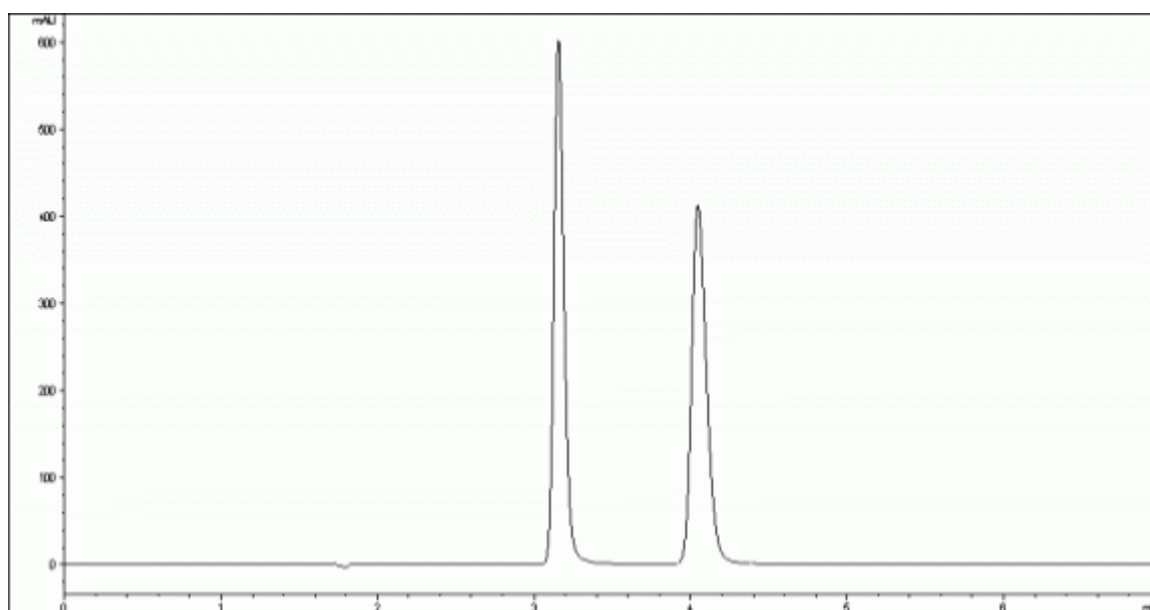

ANALYTICAL REPORT

| | | | |
|-----------------------------|---|--|--|
| Column Trade Name | CHIRALPAK® AY-3 | | |
| Column Size (Particle Size) | 4.6 x 150 mm (3 µm) | | |
| Mobile Phase Composition | n-hexane / ethanol / diethylamine = 80 / 20 / 0.1 | | |
| Chromatographic Mode | Normal Phase | | |

| | | | |
|-----------------------|---------------|----------------------|----|
| Flow rate (mL/min) | 1 | Temperature (°C) | 25 |
| Injection Volume (µL) | 0.005 (mg) | Sample Concentration | |
| Detection | UV-VIS 270 nm | | |
| Elution Order | | | |

| No. | Rt | K' | Alpha Value | Resolution |
|-----|-----|------|-------------|------------|
| 1 | 3.2 | 0.78 | 1 - 2 | 6.29 |
| 2 | 4.1 | 1.29 | | |