


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

6-Hydroxy-2,5,7,8-tetramethylchromane-2-carboxylic acid

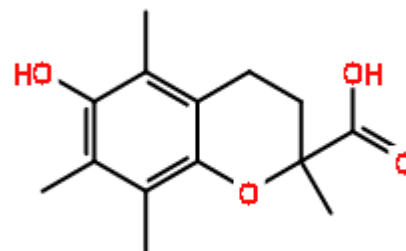

## CHROMATOGRAM

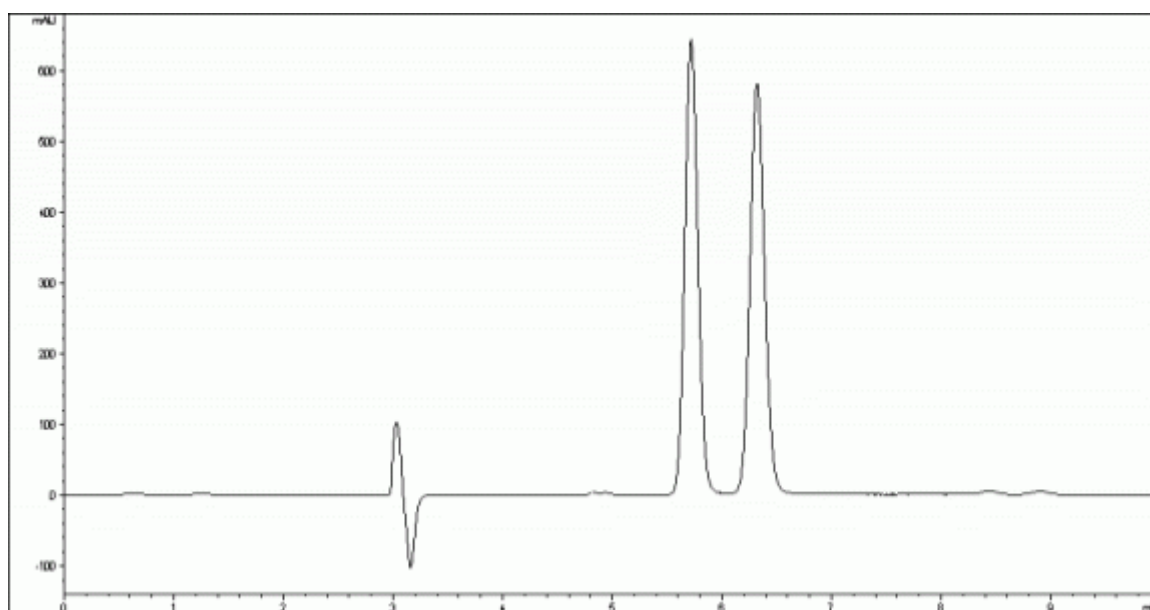

## ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® AZ-H 4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / trifluoroacetic acid = 90 / 10 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.7	0.93	1 - 2 1.22	2.60
2	6.3	1.13		