


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

1-Phenyl-2-propyn-1-ol

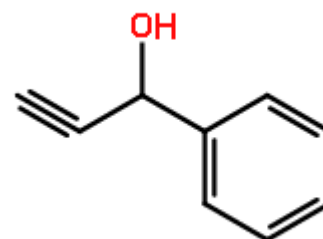

CHROMATOGRAM

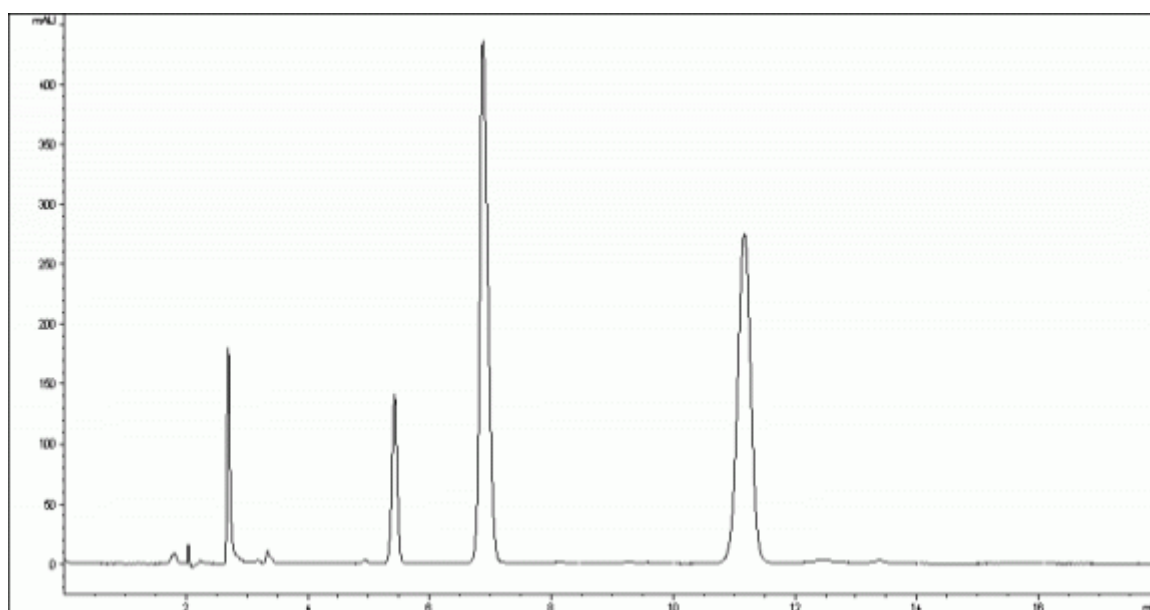

ANALYTICAL REPORT

Column Trade Name	CHIRALCEL® OJ-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	n-hexane / 2-propanol = 80 / 20		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.9	2.87	1 - 2	13.56
2	11.2	5.29		