


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Metalaxyl

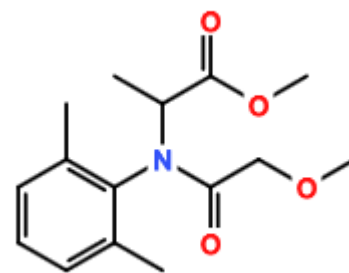

CHROMATOGRAM

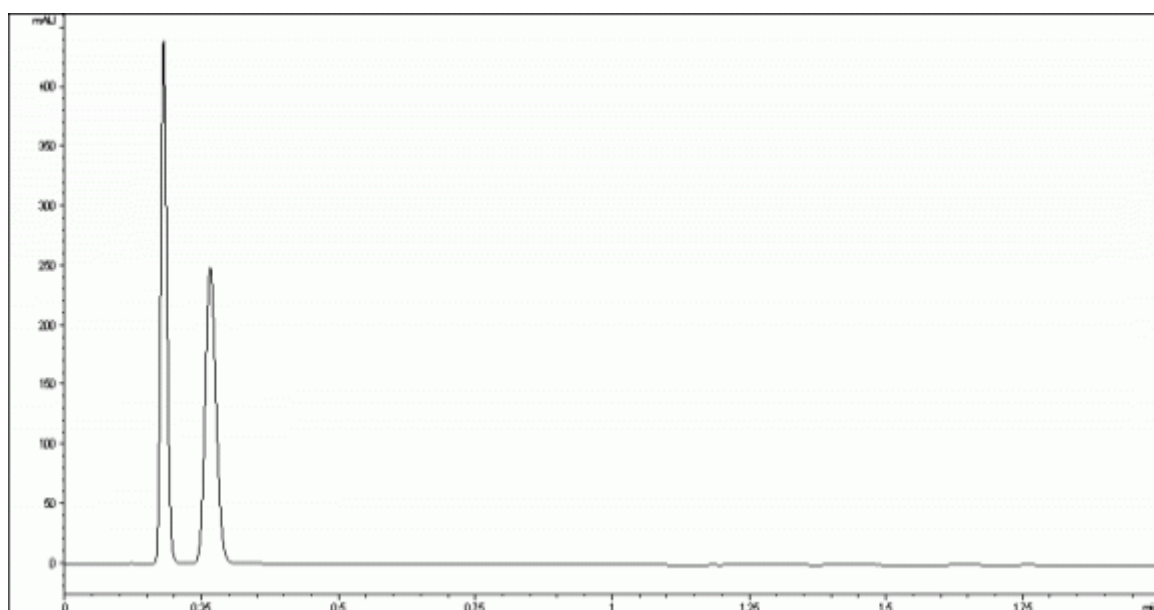

ANALYTICAL REPORT

Column Trade Name	CHIRALCEL® OD-3		
Column Size (Particle Size)	4.6 x 50 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol = 70 / 30		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	5	Temperature (°C)	40
Injection Volume (µL)	0.001 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	0.2	0.62	1 - 2	2.23
2	0.3	1.38		
				3.00