


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

1-Indanol

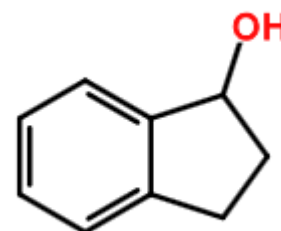

## CHROMATOGRAM

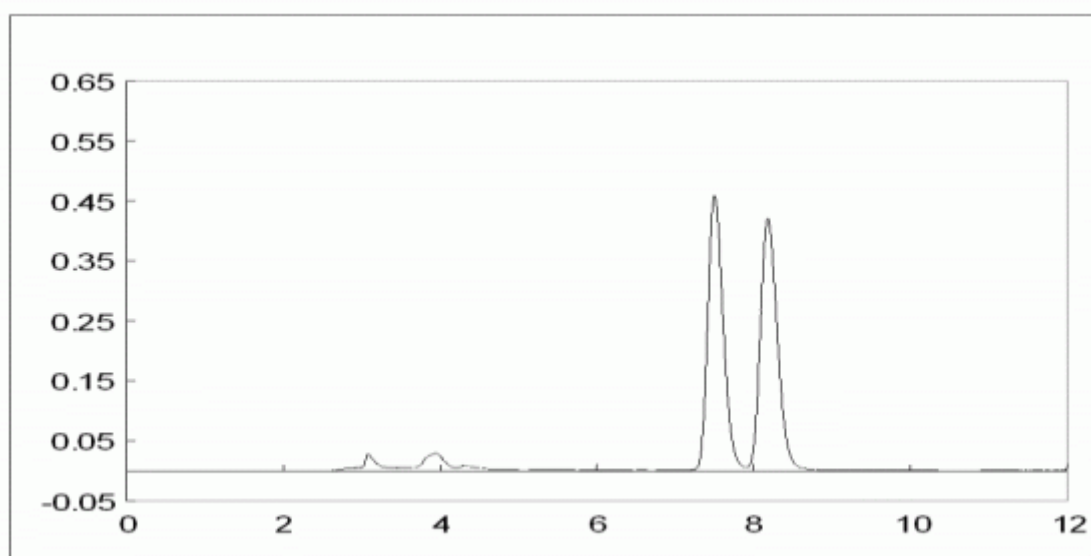

## ANALYTICAL REPORT

Column Trade Name	CHIRALCEL® OD-H		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-heptane / 2-propanol = 95 / 5		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.5	1.49	1 - 2	1.80
2	8.5	1.73		