


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Labetalol

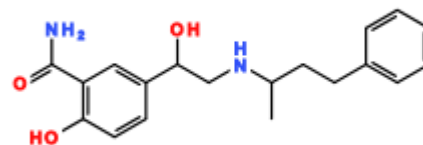

## CHROMATOGRAM

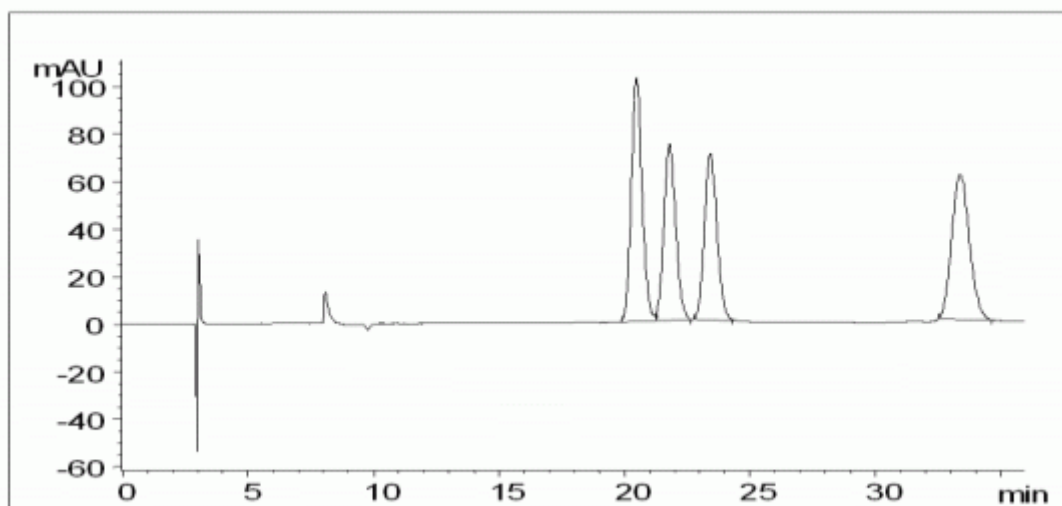

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AD-H		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / ethanolamine = 90 / 10 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	20.5	5.80	1 - 2	1.60
2	21.8	6.30	2 - 3	1.90
3	23.4	6.80	3 - 4	9.10
4	33.3	10.10		