


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Mandelic amide

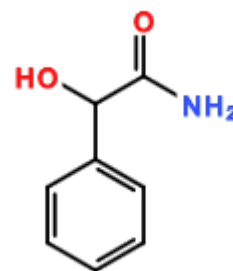

## CHROMATOGRAM

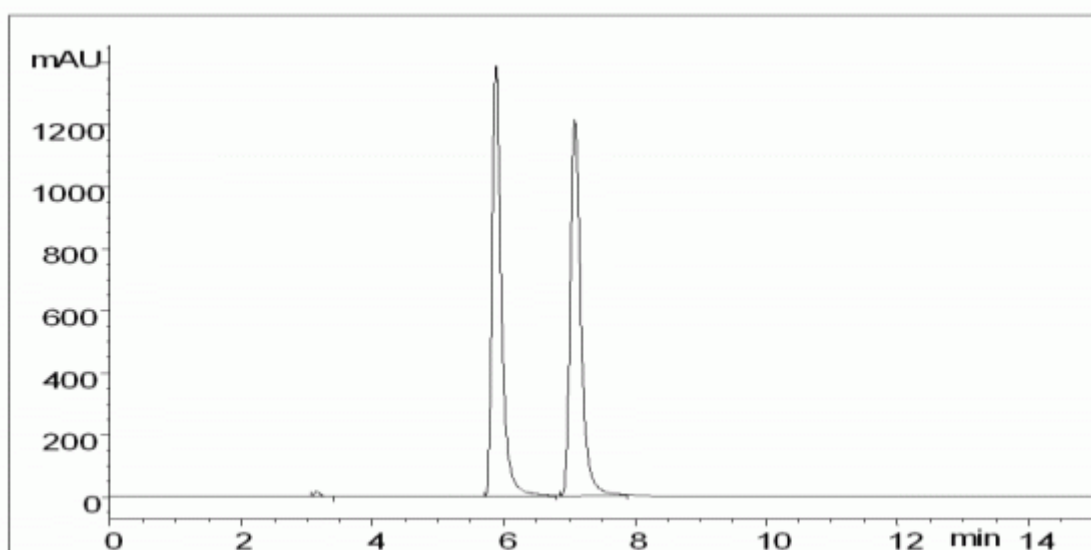

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AD-H		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-heptane / 2-propanol = 80 / 20		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.9	0.96	1 - 2	4.50
2	7.1	1.36		