


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Nicardipine

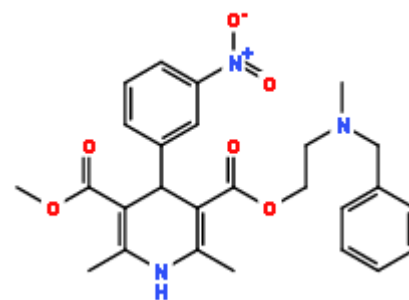

CHROMATOGRAM

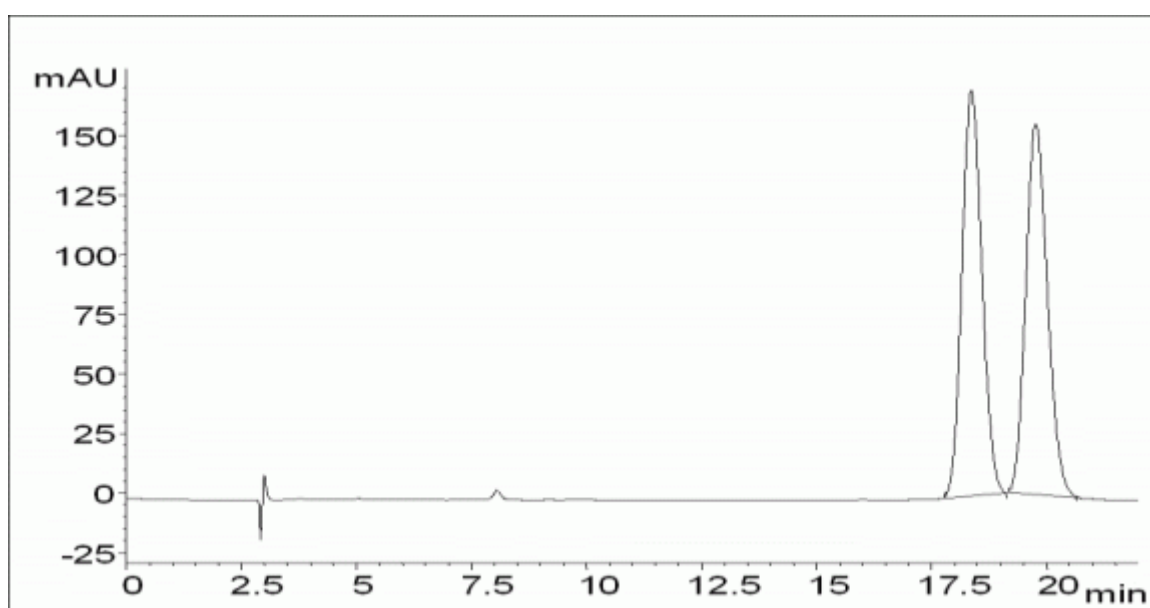

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AD-H		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / 2-propanol / diethylamine = 93 / 7 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 275 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	18.4	5.12	1 - 2	1.70
2	19.8	5.58		