


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Ofloxacin

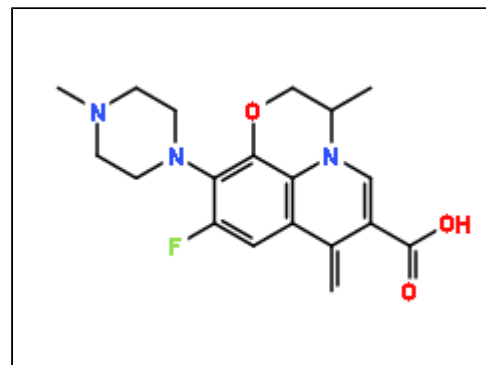

CHROMATOGRAM

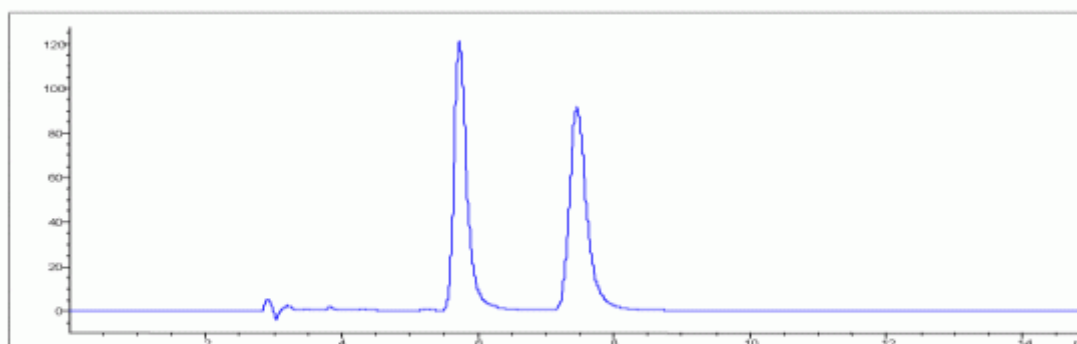

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AS-H		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / acetic acid / ethylenediamine = 50 / 50 / 0.3 / 0.3		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	40
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 300 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.7	0.91	1 - 2	4.54
2	7.5	1.51		