


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

1-Phenylethanol

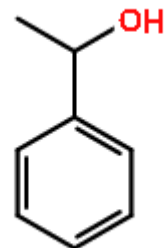

CHROMATOGRAM

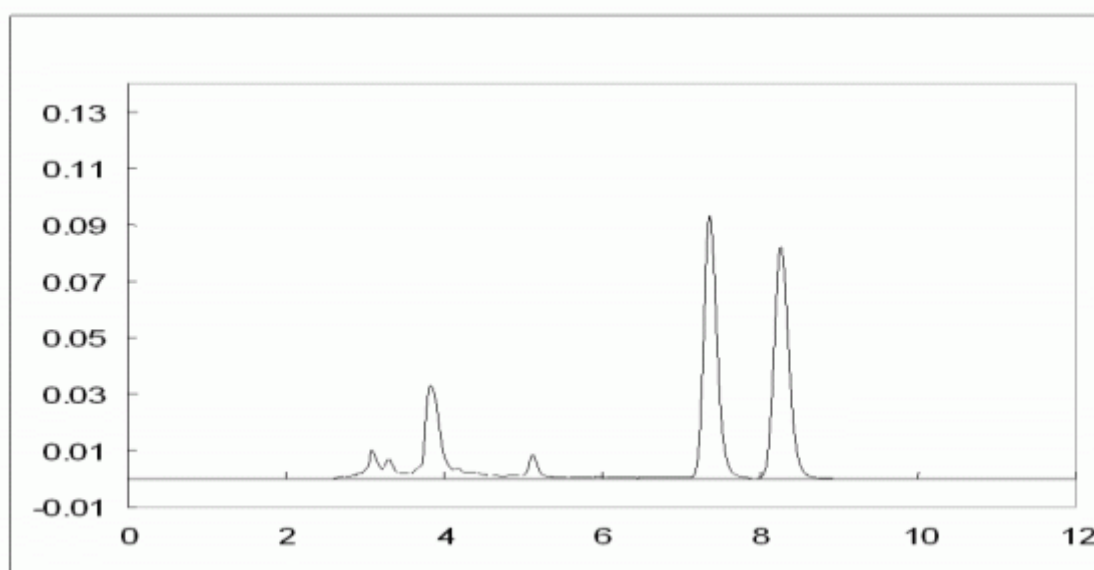

ANALYTICAL REPORT

Column Trade Name	CHIRALCEL [®] OD-H		
Column Size (Particle Size)	4.6 x 250 mm (5 μ m)		
Mobile Phase Composition	n-heptane / 2-propanol = 95 / 5		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature ($^{\circ}$ C)	25
Injection Volume (μ L)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.3	1.43	1 - 2	2.80
2	8.2	1.74		