


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

trans-Stilbene Oxide

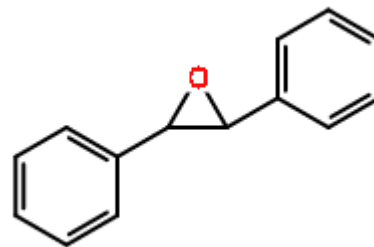

## CHROMATOGRAM

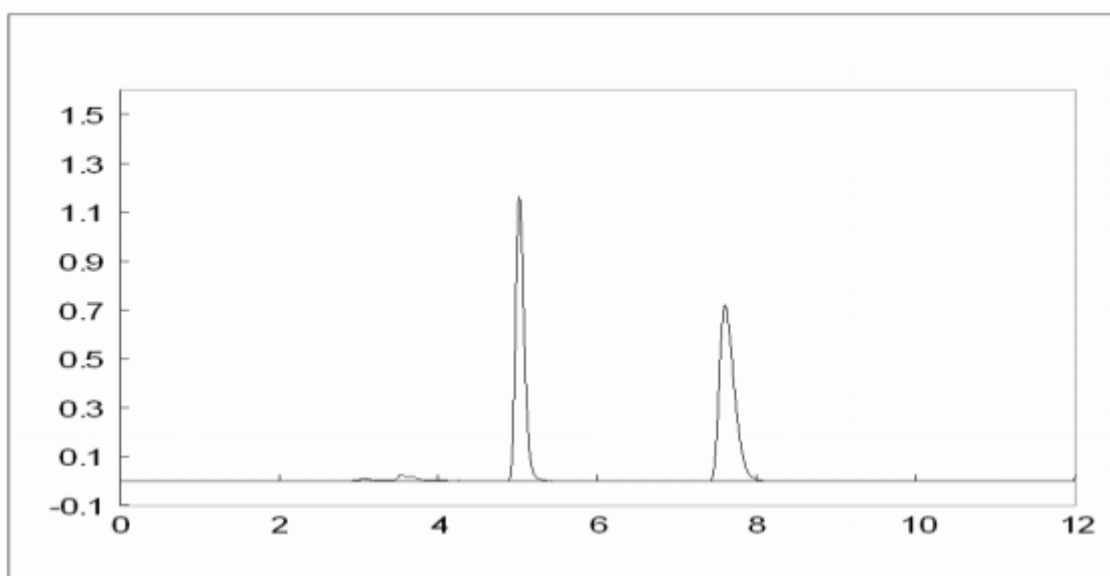

## ANALYTICAL REPORT

Column Trade Name	CHIRALCEL <sup>®</sup> OD-H		
Column Size (Particle Size)	4.6 x 250 mm (5 $\mu$ m)		
Mobile Phase Composition	n-heptane / 2-propanol = 90 / 10		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature ( $^{\circ}$ C)	25
Injection Volume ( $\mu$ L)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.0	0.65	1 - 2	9.52
2	7.6	1.50		