


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

(±)-N,N'-Bis(alpha-methylbenzyl) sulfamide

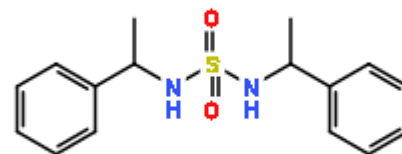

CHROMATOGRAM

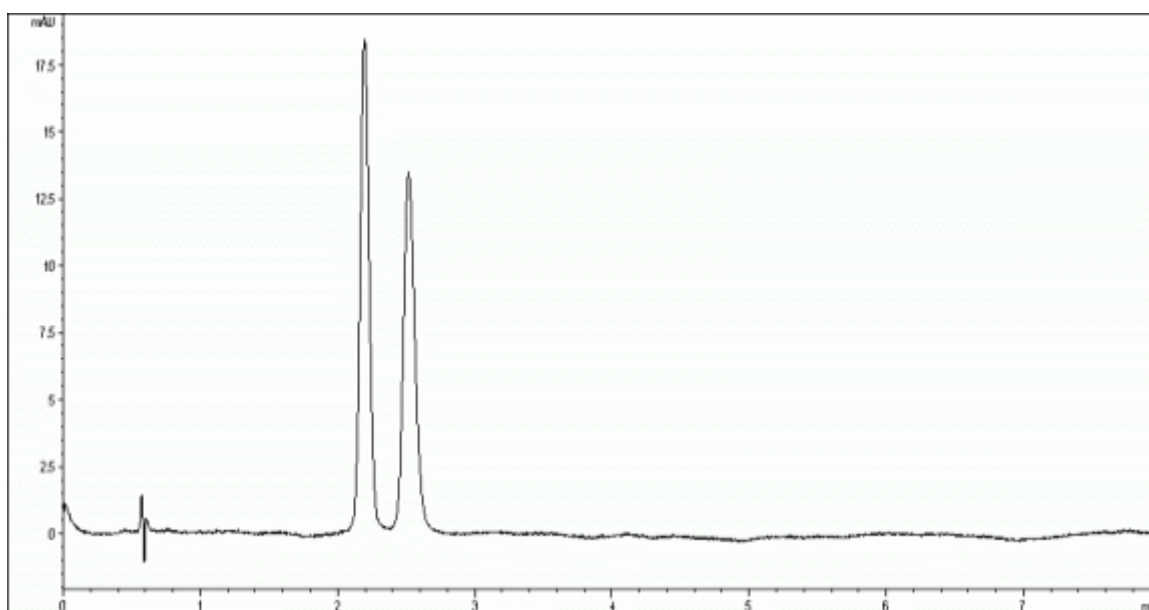

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 50 mm (3 µm)		
Mobile Phase Composition	water / acetonitrile = 50 / 50		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	2.2	2.68	1 - 2	2.39
2	2.5	3.22		