


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Cloprostenol

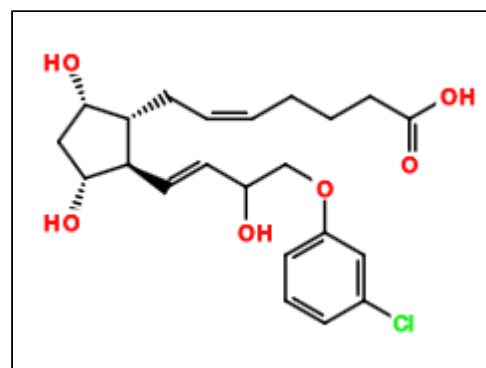

## CHROMATOGRAM

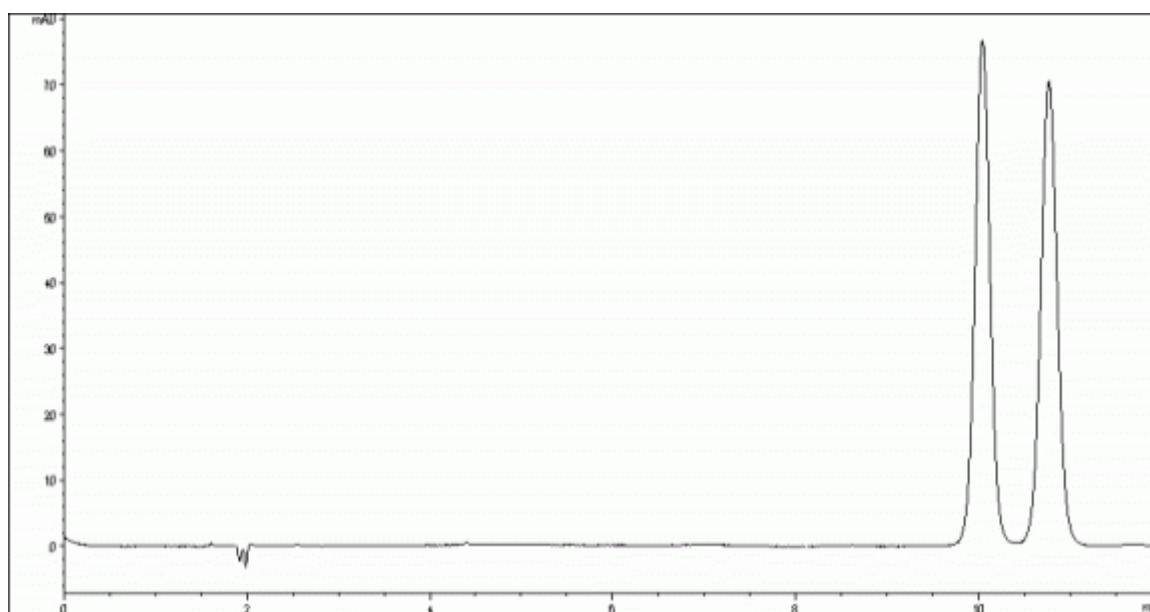

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	trifluoroacetic acid aq. pH2.0 / acetonitrile = 70 / 30		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	10.1	4.74	1 - 2	1.09
2	10.8	5.16		
				2.33